
## Bian Stone far infrared warm compress physiotherapy Instruction


### Product introduction

Dr.Rock Air is a Bian Stone heating shock body lightweight massage gun. The Bian Stone warm head can emit far infrared and ultrasonic pulses that are beneficial to the human body. Significant results have been achieved for the relief of sciatica, periarthritis, muscle soreness, and meridian dredging throughout the body. It can be changed by different Bian Stone massage heads.

### Schematic diagram of product structure
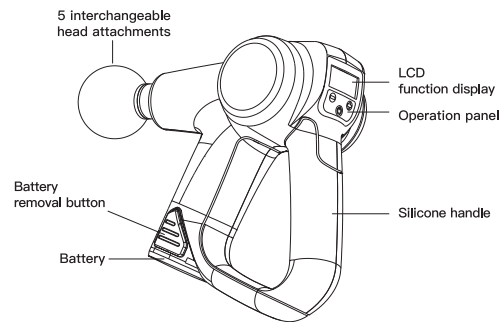
### Installation method of battery and physiotherapy head

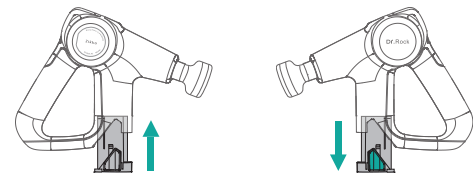
#### Battery installation method

★The convex groove of the main machine is aligned with electricity The groove of the pool is pushed into place by force.


#### Battery disassembly method

★Press and hold the battery removal key. Pull out the battery with vertical force.

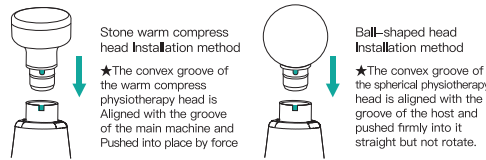**Note:** 1. Bian Stone warming physiotherapy head and other physiotherapy heads have installation limits, which are inserted after alignment; 2. Please note that the battery cannot be rotated when dismounting. Press the dismounting button and insert it vertically.

### Instructions for use

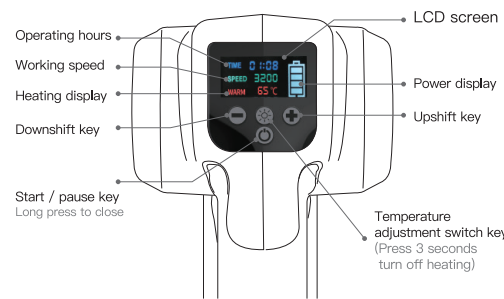- After Dr.Rock Air inserted the battery, short press the "⏻" key, the device is on warm compress;
- When the device is on, the LCD display shows the current Dr.Rock Air startup status:
  - "Battery graph" shows the current battery level;
  - "TIME 88:88" shows working time;
  - "SPEED 8888" shows the working speed;
  - "WARM 88°C/0FF" shows that the temperature is turned on (40/50/65) or closed.
- After Dr.Rock Air is started, short press the "⏻" key momentarily to pause or resume work, long press for turn off
- When Dr. Rock Air start, you can press the "⏻" key to start the warm compress function. The temperature has 3 levels, the initial low temperature is 40°C, the second comfortable temperature is 50°C, and the third gear is 65°C. The proper temperature can be adjusted reasonably according to one's own physique and preferences. Press and hold the "⏻" button to turn off the warm compress function. (Note: When using warm compress physiotherapy, it is necessary to customize the matching Bianstone physiotherapy head with the instrument);
- When Dr. Rock Air working, short press "⏻" / "⏻" key to adjust the working mode according to your own physique and preferences;
- Dr. Rock Air has intelligent work protection. When you stop use about 3 minutes and it is not in warm compress mode, it will automatically turn off, enter the sleep status. If you need to use it, you need to restart it. (Note: When the warm compress mode is on, it will not enter sleep status.)

### Characteristics and Efficacy of Physiotherapy Head

<p><b>Bian Stone flat warming head</b></p> <p>40~65°C Specialized treatment to relieve sciatica, frozen shoulder, Joint pain, muscle aches, swelling, Systemic body meridians dredge and other pain.</p>	<p><b>Bian Stone cone warming head</b></p> <p>40~65°C It is suitable for impacting deep tissues, such as meridians, joints, palms and soles.</p>
<p><b>Bian Stone ball warming head</b></p> <p>40~65°C It is suitable for massaging large muscle groups such as arms, waist, back, buttocks, thigh, calf, and etc.</p>	<p><b>Ball-shaped head</b></p> <p>It is suitable for massaging various muscle parts and other large muscle groups.</p> <p><b>U-shaped head</b></p> <p>It is suitable for massaging neck, spine and Achilles tendon.</p>

Accelerating warm-up effect  
Maintain muscle activity  
Fast decomposition of lactic acid  
Deep fatigue  
Relieve muscle tension  
Rapidly release lactic acid  
Vibration relaxation of muscle fascia  
Stretch muscle lines  
Improve muscle balance

Activating muscle before exercise state during exercise  
elimination of fatigue after exercise  
effective muscle relaxation pain and stiffness  
accelerate muscle recovery  
relaxation of muscle fatigue aid in fitness plasticity  
assist all kinds of training

## Bian Stone far infrared warm compress physiotherapy

Warm compress therapy is an important method in stone therapy. Its origin is to use baked hot stone to rest on the affected part of the disease or massage, which can relieve pain and treat the disease. This therapy is called stone warm compress. Electrothermal stone therapy is a kind of warm compress, i.e. current is used to heat the stone through a special instrument. The far infrared electromagnetic wave radiation heat energy generated by the heated stone can be conducted in channels and collaterals to achieve the purpose of treating various diseases. It is a clever combination of traditional stone therapy and modern high technology. At the same time, stone has a strange energy field, which acts on the surface of human skin to generate extremely far infrared radiation with extremely wide frequency band and far infrared frequency ranging from 7 to 20 microns. Every time stone rubs the human body, it can generate ultrasonic pulses beneficial to health. The average number of ultrasonic pulses can reach 3708, with the frequency range of 2-2 million Hz. Therefore, stone is the first of all materials.

Dr.Rock series Bian Stone far infrared warm compress physiotherapy has similar effect to that of warm compress therapy. It is a compound therapy of stone hyperthermia and massage, and has excellent curative effects of promoting blood circulation, removing blood stasis and dredging channels and collaterals. After thousands of years of practice, warm compress therapy has proved to have a wide range of health care functions. Clinical evidence shows that electrothermal stone therapy has good effects of warming channels and dredging veins, eliminating arthralgia, relieving spasm and sedation, and can be used for various pains, relieving muscle spasm, digestive system diseases, and etc. It also shows great advantages in beauty treatment, weight loss, insomnia adjustment and fatigue elimination.

### Massage impact therapy

#### Impact on skin and subcutaneous tissue

Loosening thickened connective tissue and fascia, breaking the barrier of releasing fascia contraction, removing adhesion and internal scar tissue after injury or surgery, increasing venous return and lymphatic circulation, and improving muscle extensibility. Especially when muscle spasm, stiffness and restriction occur, the reflex contraction of muscle fibers is punctured, and the range of motion is increased by relieving stiffness of joints.

#### Influence on central and sympathetic nervous system

Stimulate local muscle receptors to relieve pain and help relax, without causing pain or inflammation.

#### Influence on muscle, preserved muscle, arthralgia

The most common causes of muscle pain are trauma, overuse, tension or tension. After injury or trauma occurs, impact therapy can accelerate the healing of inflammatory regions. Similarly, when some parts of the body are overused or under tension, impact therapy can accelerate blood circulation in affected parts and reduce potential tension.

#### Exercise rehabilitation

Relaxation physiotherapy is helpful to recover the body from the pressure of strenuous exercise and to adjust and rebuild training. Rehabilitation will remove metabolic waste and better nutrient cells, normalize and enlarge tissue elasticity, treat trigger points and accelerate wound healing.

#### Recovery after operation and injury

After trauma or surgical repair, the stimulation of impact therapy can eliminate faults and avoid formation and adhesion of scar tissue during the repair process of inflammatory activity of the body.

#### Relax

Impact physiotherapy can repair the body, increase the patient's body comfort and improve the overall well-being.

## Indications and contraindications

Dr.Rock series deep penetration can relieve pain, bursitis, muscle spasm, neuralgia and other musculoskeletal diseases caused by joint sprain in a short time. Almost all soft tissue diseases are related to trauma, overuse of muscles or abnormal posture. Dr.Rock can be used to treat scar tissue involved in most of these diseases. When using Dr.Rock, excessive pressure and long-term physical therapy should be avoided. Sensitive areas such as head, face and cervical vertebra should also be avoided. Physical therapy should be stopped immediately if there is inflammation, swelling or pain aggravation. The current research is about the impact of vibration training on specific medical diseases. This is likely to reduce the types of contraindications mentioned above. Practical experience shows that in many cases, the combination of vibration training with physiotherapy program is also a physiotherapy suggestion. This must be done on the advice or accompanied by a doctor, expert or physiotherapist.

### Indication

- ✓ Pain and spasm due to muscle contusion, sprain and strain
- ✓ It can help edema fluid flow in swollen parts,
- ✓ Relax thickened connective tissue and fascia,
- ✓ Reduce lactic acid accumulation in muscles, increase joint activity,
- ✓ Eliminate muscle fatigue and
- ✓ Sciatica (Bian stone hot compress physiotherapy) and
- ✓ Scapulohumeral periarthritis (Bian stone hot compress physiotherapy)
- ✓ Periarthritis (Bian stone hot compress physiotherapy)

### Contraindications

- Aneurysm, hemorrhage, heart disease with blood thinner;
- cardiac pacemaker or defibrillator;
- Pregnancy, cancer;
- It cannot be used within 3 inches of internal fixation when joint internal fixation less than 90 days,
- Sensitive areas: head, face, cervical spine, Spinal vertebrae
- Muscles cover less locations near bones:  
Tibia, back of foot, back of hand

## Matters needing attention in physical therapy

### After injury

- The first physical therapy after injury should be based on the patient's satisfaction and comfort level
- Rom evaluation
- Determine the position and surrounding muscles according to the physical therapy plan
- Start the instrument and set it at medium and low strength.
- The applied pressure is the weight of the instrument itself

### After surgery

- The first physical therapy after operation gradually increases the time and intensity of physical therapy according to the patient's reaction and feeling.
- Determine physiotherapy site

- Ask if the patient has open wound, suture needle, drainage and internal fixation. (If yes, do not stimulate these parts)
- Start the instrument and set it at the lowest gear.
- Cover the physiotherapy area with towel.
- Apply pressure to the physiotherapy site with appropriate amount of instruments.
- Physical therapy is limited to the area around the surgical site.

### Lymphatic drainage

- Minor and systemic physical therapy is required.
- During physiotherapy, the limb moves slowly from the distal end to the proximal end, i.e. towards the direction of cardiac vein reflux

### effect

- Reduce swelling and pain
- Eliminate scar tissue after operation and injury
- Restore the elasticity and extensibility of muscles
- Trigger therapy
- Promote blood circulation and lymph reflux
- Myofascial relaxation

## Maintenance and cleaning

### Maintenance, cleaning and storage:

- Clean the equipment with a slightly wetter towel and dry it with a soft cloth.
- Before and after use, please wipe and disinfect the stone physiotherapy head with medical alcohol.
- After each use of the product, turn off the main power switch.
- When storing or travelling, remove the battery and put it in the paper box provided.

### Charging:

- The charger standard for this product must be used for charging, otherwise our company will not be responsible for any consequences arising therefrom.
- Before the first use, the battery needs to be fully charged before use
- When the battery is charging, the indicator light on the charger is red, and the indicator light is green when fully charged

### Warnings and reminders:

- Do not use this machine when charging
- Do not charge the product for a long time, such as more than 12 hours.
- Do not throw, expose or drench the product violently. Keep away from liquid or heat source and keep the vent free of ash, dust and debris
- Do not remove screws or attempt to disassemble products.
- After each use of the product, turn off the main power switch.
- For adults only. Please consult your doctor before using if you are injured.

## Safety precautions

- For adults only;
- It can only be used on dry and clean skin surface of body parts, and cannot be used through clothes. It is light. Gently press and move on the skin, each part for about 60 seconds;
- Use only on soft tissues of the body without causing pain or discomfort. Not on the head or Use on any hard, bony part of the body;
- Please use the physiotherapy head that can produce the most ideal effect.
- Both controlled depression and pressure may cause bruises. Check physical therapy area frequently. Stop use immediately if there are signs of pain or discomfort;

- Do not let fingers, hair or other parts of the body contact the shaft and back of the vibrating head to avoid being caught.
- Do not drop Dr.Rock or use it improperly.
- For the first time, please use it after full charge. Please use the standard original charger to charge Dr.Rock, do not charge for a long time, just fill it up.
- Carefully check the equipment and battery before each use.
- Please do not change or modify Dr.Rock's operation in any way.
- Never allow Dr.Rock to run or recharge in an unattended state.

## Packaging accessories

- Carrying case x1
- Manual x1
- Battery x2
- Charger x1
- UK plug x1
- Dr.Rock Air x1
- Use tutorial x1
- Quick guide x1
- CN/US plug x1
- EU plug x1
- Bian Stone flat head x1
- Bian Stone cone head x1
- Ball-shaped head x1
- U-shaped head x1

## Specifications and parameters

Model: H-MG300  
Name: Zikko Dr.Rock Air Bian stone far infrared warm compress physiotherapy  
Massage mode: hot compress, uniform speed  
Gear / speed: 900 ~ 3200 rpm  
Speed gear: up / down, 200 rpm / click to increase or decrease  
Power: ≤60W  
Bian stone heating temperature: first gear 40 °C, second gear 50 °C, third gear 65 °C (heat compress heat temperature)  
Charger input voltage: 110 ~ 240V / 50 ~ 60Hz  
Charger output voltage: DC 12V / 1A  
Battery: 12V / 2600mAh, 31.2Wh, lithium battery  
Duration: up to about 6 hours x2  
Charging time: about 150 minutes  
Size: 226\*174\*104mm Weight: 900 g Color: Gray

## Recommendations for handling abnormal situations

If the following abnormal situations occur when using this product, please refer to the solution:

Abnormal situation	Abnormal cause	Exception handling measures
The function cannot be started	The battery may be low on power.	Please charge
Automatic shutdown in use	The battery power may be insufficient	Please charge
warming physiotherapy head not hot enough after opening	Physiotherapy head warms up in two minutes	normal
warming physiotherapy head Temperature is too high	After warm mode is turned on Equilibrium temperature with human body	Turn off the warm compress mode for 3 minutes and then turn it back on for normal use
After 3 minutes of instrument pause (With warm pack off) The display goes off	Switch to smart power saving mode	Short press the start key to return to the previous working mode

Note: If there are other abnormal conditions, please contact customer service in time.

## After-sales service

### One-year warranty

The quality of the whole machine is guaranteed for one year from the date of purchase, and the accessories attached are guaranteed for six months. Due to the quality problems caused by the manufacturing defects of the products, they cannot be used normally, and our company will be responsible for free maintenance, replacement or replacement of necessary components.

- Damage caused by accidents, abuse, improper handling or transportation;
- Repair products without permission;
- The product is not maintained or used according to the instructions.
- Loss is higher than product cost;
- Delivery quality problems due to abnormal storage or protection conditions in the customer's workplace;
- Unable to provide valid purchase voucher;
- Products purchased through unofficial or authorized channels

Some states and countries do not allow restrictions on fault compensation, so the above restrictions may not be available to you. This warranty covers specific legal rights, but other rights may vary by country and state. This warranty service is provided by Zikko (Hong Kong) Digital Technology Co., Ltd. and Zikko (China) Co., Ltd. If you are a customer, please contact the physical store, distributor or retailer where you purchased the product directly to learn about the warranty and return procedures.

## Register/註冊

In order to guarantee your benefits, please register by zikko official website www.zikkotech.com within 90 days from the date of purchase, to gain more professional services and supports from Zikko. pls refer the detail customer service description on our website.

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#### Register/注册

In order to guarantee your benefits, please register by zikko official website www.zikkotech.com within 90 days from the date of purchase, to gain more professional services and supports from Zikko. (to refer the detail customer service description on our website. 為保障購買Zikko (即刻) 用戶的權益，請用戶於購買之日起90天內在Zikko (即刻) 官網 www.zikkotech.com 進行產品註冊，以獲得Zikko (即刻) 更多專業的服務與支持！具體售後細節詳見網頁。

Scan a QR Code 掃碼獲取更多資訊  
www.zikkotech.com

8