**Dr.Rock Pad**  
BianStone Far Infrared Heating Therapy SitPad  
砵石遠紅外線熱護坐墊

**User's manual**  
使用說明書

Scan QR Code  
請掃描QR碼

**Zikko**  
Room 210, 2/F, Block A, Kai Fu Chase 3, Sheung Shui  
Sheung Shui, N.T., Hong Kong

**CE** **FC**

Designed by Zikko in Hong Kong  
www.zikkotech.com

Copyright © 2022 Zikko (HK) digital technology Co., Ltd. All rights reserved.  
Zikko & Dr.Rock are registered trademarks of Zikko (HK) digital technology Co., Ltd. All other trademarks are the property of their respective owners.

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Sign/Seal Here :  
經銷商簽章:

**Dr.Rock BianStone Far Infrared Heating Therapy SitPad**

**Product Introduction**

Dr.Rock BianStone Far Infrared Heating Therapy SitPad is made of pure natural premium BianStones of red colour. Electronic constant temperature heating of BianStones can emit far-infrared rays beneficial to the human body, accelerate and improve the blood micro-circulation of the human body, activate cells, and achieve thermal curative effect. It can help relieve body fatigue and pain, relax muscles, improve acidic physique and sleep quality and other warm care and warmth. It is a pure natural health physiotherapy method for warming the whole body. Long-term use it can promote metabolism and enhance human immunity.

**Eight characteristics**

Natural premium BianStone

Graphene heating film

Far-infrared wavelength 6-15 microns

Anti-burning Anti-leakage

Full touch designed

Heating 20~60°C

Slim and easy to be rolled


Soft leather

**Product Structure**

- Nylon edging
- High-density mesh cloth
- Natural BianStone
- Bamboo charcoal fiber cotton layer
- Graphene heating film
- Gathering energy and slow heating layer
- Soft leather


**Product usage diagram**

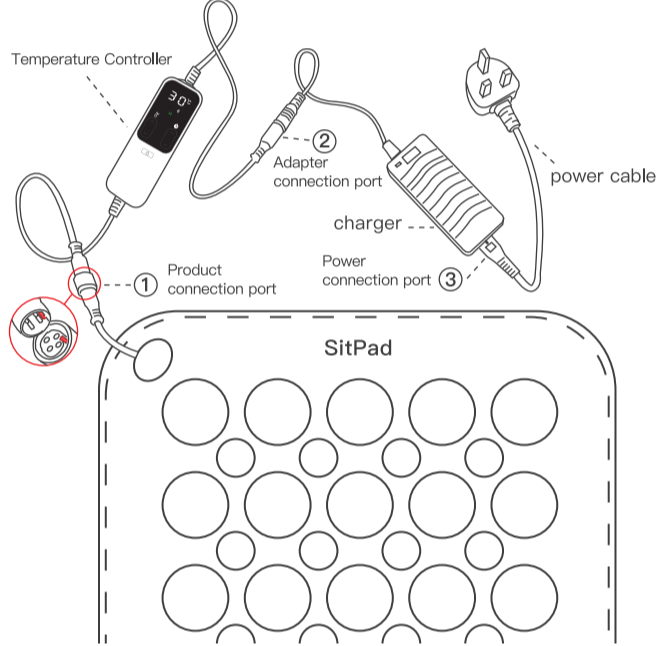**Product connection diagram:**

- Connect the four-pin connector (male) with the thermostat's four-pin connector (female) according to the limit slot and arrow to confirm the butt connection;
- Connect the DC connector (female) at the other end of the thermostat and the DC connector (male) of the power supply end to the limit slot attached to it and rotate and fix it;
- Insert the end of the matching power cord into the power connector.

**Controller use and troubleshooting**

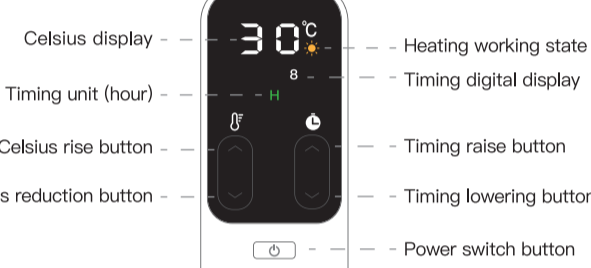**mode**

mode	set up	Automatic heating time
A.Warm therapy 1	20~45°C	Can be set manually for 2~12 hours
B.Warm therapy 2	46~50°C	Heating for 4 hours, automatically return to 37°C
C.Hot therapy	51~60°C	Heating for 4 hours, automatically return to 37°C

**6. Thermostat display troubleshooting**

(display)	cause of issue	Troubleshooting
E1	The thermostat is not connected to the product or the connection is in poor contact. ② The temperature sensor probe fails.	1. Check the connection line, whether the connection is normal, and it needs to be re-wired. 2. If the product is still abnormal, please contact after-sales service.
EH	High temperature error or temperature out of protection	Turn off the power, turn on the product again after cooling to confirm, if it shows that the product is abnormal, please contact after-sales

**The product soothes and improves the following people**

- Hemorrhoids, insomnia, long-term fatigue, deficiency and cold physique;
- Female menstrual period, cold hands and feet, edema, obesity, sub-health
- Back pain, shoulder and neck pain, numbness of hands and feet;
- Sports injuries.

**Physiotherapy mode**

- Warm Therapy:** Set the temperature to 30~50°C, for daily conditioning of the body.
- It is suitable for cold hands and feet, female physical discomfort, insomnia, improvement of acidic physique, etc. It can promote blood circulation and obtain a sense of comfort. It is recommended to use it for 10 hours/each time.
- Hot Therapy:** The temperature is about 51~60°C, and it is used for hyperthermia when needed.
- Used for sports injuries, backache and bone pain, shoulder and neck pain, numbness in hands and feet. Can accelerate blood circulation, perspiration and detoxification, effectively relieve and improve pain. Recommendations for use: 2 hours/each time.
- Cold Therapy:** The power heating switch is turned off (room temperature) for natural care at normal room temperature.
- It is cool and comfortable to use in summer; long-term use can make the body absorb the far infrared rays naturally emitted by the BianStone, which is beneficial to the body, which is beneficial to maintain and improve the normal operation of the body's microcirculation system. Recommendations for use: long-term.

**Specifications**

Model : SP450  
Product name : Zikko | Dr.Rock BianStone Far Infrared Heating Therapy SitPad  
Size : 45cm x 45cm / Weight : 1.2 KG  
BianStones quantity : Large grains 5 rows x 5 grains + small grains 4 rows x 4 grains = 41 grains

Therapy mode : cold therapy (indoor temperature is not heated), warm therapy, hot therapy

Adapter input voltage : 100-240V~50/60Hz 2.0A  
Adapter output voltage/current : 24V 2.75A  
Thermal cushion power : 72W max  
SitPad Current and Voltage : Human body safety low voltage max 24V/3A  
Heating device : Graphene heating film  
Adjustable range of temperature controller : 20~60°C  
The temperature controller can set the timing : 2~12 hours  
Heating temperature : 20~60°C up and down button adjustment  
Accessories : soft leather, fireproof mesh, double-layer felt cloth

**Package Contents**

\*BianStone SitPad x1 \*Manual book x1 Power adapter x1  
\*Temperature controller x1 \*Product sheath x1

**Warranty description**

Zikko provides quality warranty services for hardware products and accessories in the packaging and defects in materials and workmanship. The warranty period of the Dr.Rock SitPad (including temperature controller) is one year from the date of purchase of the product, and the included storage bag is guaranteed for 6 months. Damages caused by normal wear or aging, accidents, improper use, alterations, etc. are not covered by this warranty. How to obtain after-sales service: Please call the Zikko authorized brand or go to the authorized retailer's purchase place or online customer service for after-sales service, and please provide the relevant purchase credentials (invoice or purchase record).

**Product maintenance methods**

You can use a wet towel to wipe the contaminated part when it is not powered on, or use a disinfectant wet towel to disinfect it regularly, and then turn it on for heating after drying.

**Precautions**

- Read the manual carefully before use and use it in accordance with the operating instructions.
- Don't let it be used alone by people who cannot take care of themselves, infants or children who are not sensitive to heat.
- It is forbidden to use in water and humid environment.
- Do not dry clean or wash. You can wipe the contaminated part with a wet towel when the power is not turned on, and then turn on the power to heat it after drying.
- Do not insert sharp objects into this product.
- Do not place bedding and other items in a small area to cover the heating area when the product is heated to prevent local overheating and causing malfunctions.
- It is forbidden to place flammable materials such as fireworks and lighters on the surface of the product to prevent fire.
- Please confirm whether the power supply is out off before going out.
- Using the heating function, the temperature can be adjusted according to the tolerance of the skin. In order to prevent skin burns, if you feel unwell, stop using it in time.

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