Dr.Rock

# Dr.Rock Mini 2 Massage Gun 按摩儀

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
### Dr.Rock Mini2 Massage Gun

#### Product introduction

Dr.Rock Mini 2 massager Gun is a wireless, light and portable vibration relaxation device that can deeply relieve muscle stiffness and soreness and increase blood circulation, activity performance, and improve the overall health of the body's soft tissues. With three-speed massage mode, it is the best muscle relaxation tool necessary for home, travel and sports.

#### Schematic diagram of product structure

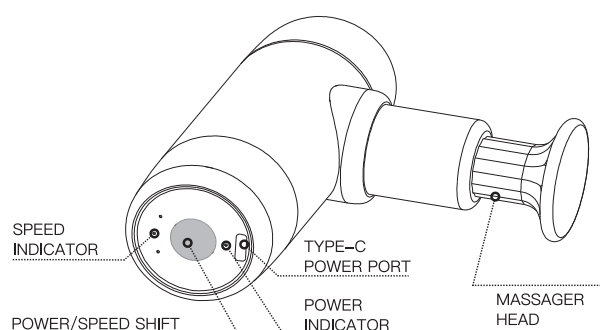
#### Massage head


#### Instructions

##### INSTALLATION STEPS


Choose the massager head to insert into the round hole
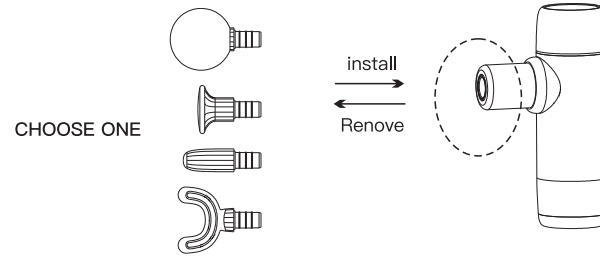
##### CHOOSE ONE
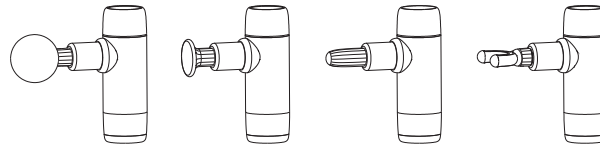
#### SWITCH
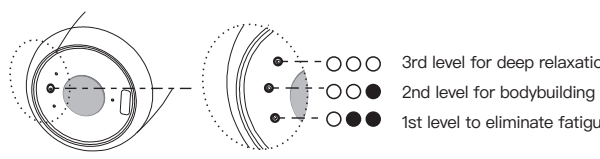
##### ON/OFF

Press and hold the power button for 2 seconds to turn on or off.


- \* It is recommended to use it after running for 5 seconds.
- \* Do not use while charging!
- \* To avoid massaging a position for a long time, each speed will automatically shut down after 15 minutes of use.

##### SPEED SHIFT

After turning on, short press the power button to switch speed.

- \* Cyclic switching of 1st, 2nd, and 3rd level.
- \* Because it's a metal, it is normal to get hot for a long time.

#### Charge

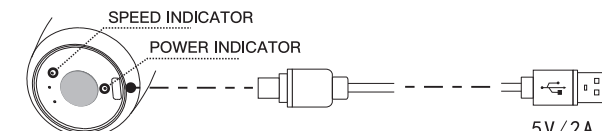POWER INDICATOR	STATE	DESCRIPTION
WHITE	FLASH	Battery 11%~30%
RED	FLASH	Battery ≤10%
RED	ALWAYS ON	Charging
WHITE	ALWAYS ON	Finished charging

SPEED INDICATOR	STATE	DESCRIPTION
WHITE	●●●	1ST LEVEL
WHITE	○○●	2ST LEVEL
WHITE	○○○	3ST LEVEL

#### Specifications

Model: H-MG110  
Name: Zikko Dr. Rock Mini2 Massage Gun  
Massage mode: uniform speed  
Speed gear: 1800/2400/3200 rpm, Click upshift  
Power: ≤30W  
Input voltage: DC 5V / 2A  
Battery: 2200mAh, 8.14Wh, lithium battery  
Duration: up to about 8-9 hours  
Charging time: about 3 hours  
Size: 84\*141\*45mm (Not including massage head)  
Weight: 460g

#### What's Included

- Dr.Rock Mini 2 x1
- User's manual x1
- Ball Massage head x1
- Cone Massage head x1
- Flat Massage head x1
- Charging Cable x1

#### NOTE

To reduce the risk of malfunction, electric shock, injury, death and equipment damage, please read this manual carefully before using this product.

- It is forbidden to use on non-dry, non-clean skin surface, you can use it across clothes, gently press and move on the skin;
- It is forbidden to use on areas above the neck and bones, and only use this product on the soft tissues of the body without pain or discomfort;
- It is forbidden to use in places where bruise or discomfort occurs;
- It is forbidden to put any object in the hot place of the product;
- It is forbidden to immerse the product in water, and do not allow water to enter the product;
- It is forbidden to put the product into the fire or near the fire source;
- It is forbidden to drop the product intentionally;
- It is forbidden to disassemble this product in any way;
- It is forbidden to insert fingers, hair or other parts of the body into the mounting hole of the massage head, it may be caught;
- Teenagers (12-18 years old) need to be accompanied by their parents. They are forbidden to use alone and children under 12 years old;
- If you need to use this product for muscle damage, please consult a doctor or use it with a professional.

#### NOTES

Without the consent of the doctor, please use this product with caution under any of the following circumstances. Pregnancy, diabetic complications (such as neuropathy or retinal damage), wearing a pacemaker, epilepsy or migraine headaches, disc herniation, spondylolisthesis, vertebrae loss or spine joint rigidity, recent surgery, etc. Anything related to your health Questions, these do not mean that you can not use this product, it is recommended to consult a doctor before use. Practice has shown that in many cases, the combination of thumping training into the treatment plan is also a treatment method, but it must be carried out with the advice or accompanying of a doctor, expert or physiotherapist.

#### After-sales service

##### One-year warranty

The quality of the whole machine is guaranteed for one year from the date of purchase, and the accessories attached are guaranteed for six months. Due to the quality problems caused by the manufacturing defects of the products, they cannot be used normally, and our company will be responsible for free maintenance, replacement or replacement of necessary components.

- Damage caused by accidents, abuse, improper handling or transportation;
- Repair products without permission;
- The product is not maintained or used according to the instructions.
- Loss is higher than product cost.
- Delivery quality problems due to abnormal storage or protection conditions in the customer's workplace;
- Unable to provide valid purchase voucher;
- Products purchased through unofficial or authorized channels

Some states and countries do not allow restrictions on fault compensation, so the above restrictions may not be available to you.

This warranty covers specific legal rights, but other rights may vary by country and state. This warranty service is provided by Zikko (Hong Kong) Digital Technology Co., Ltd. and Zikko (China) Co., Ltd. If you are a customer, please contact the physical store, distributor or retailer where you purchased the product directly to learn about the warranty and return procedures.