## Shenzhen Breo Technology Co., Ltd


19/F, Hisense Southern Building, No.1777 Chuangye Road,  
Nanshan District, Shenzhen,  
Guangdong, CN, 518056

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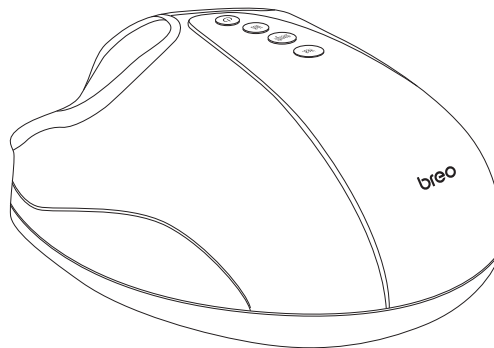**FM D191**

## **Introduction**

Thank you for purchasing Breo Foot Massager. Please read this manual carefully to ensure safe and proper use. If you need any assistance, please contact our support team with the product model number and order number.

## **Inside the package:**

- 1 x Foot Massager
- 1 x User Manual
- 1 x Adapter
- 2 x Replaceable Liners

Indications for use:

Foot Massager is intended to relieve minor foot muscle aches and pains.

Product Specifications

Model Number	FM D191
Name	Foot Massager
Power Supply	Input:100-240VAC 50/60Hz 0.8A; Output:12.0V== 2.5A
Rated Power	30W
Rated Current	2.5A
Default Time	15 min
N.W.	3500g
Dimensions	L418×W359×H 230 (mm)

Contraindications

Do not use this product if:

1. You are pregnant.
2. Your feet are injured.
3. You have a pacemaker or other medical implants.
4. You have skin infections or any other skin diseases.
5. You suffered or have been suffering from heart disease.

Precautions

1. This device does not substitute health care professionals' advice and treatment. It is not for diagnosis, therapy or other medical purposes.
2. People with impaired motor skills or suffer from physical and/or mental disorder should use this product under close supervision.
3. If at any time soreness, rash, or abnormal pain occurs, discontinue use and consult your physician.
4. The foot size at maximum for the product: 44 (EU size).
5. You are recommended to use the product for 2-3 times per day.
6. The temperature of the heating function at maximum: 55°C / 131°F.
7. The appliance has a heated surface. Persons insensitive to heat must be careful when using the appliance.
8. Consult a medical professional before use if you have any concerns regarding your health.


Troubleshooting:

Problem	Suggested Solution
The product doesn't work.	Restart the product after it cools down.
The massage intensity is too strong or too weak while using.	Adjust the angle until the intensity is moderate.
	Do not stand on or in appliance. Use only while seated.
	Put on socks before use if the intensity is too strong for you.


Any problem ,don' t hesitate to contact us for help:

Email: support@breocare.com 1-866-546-3587(United States)

Product Diagram

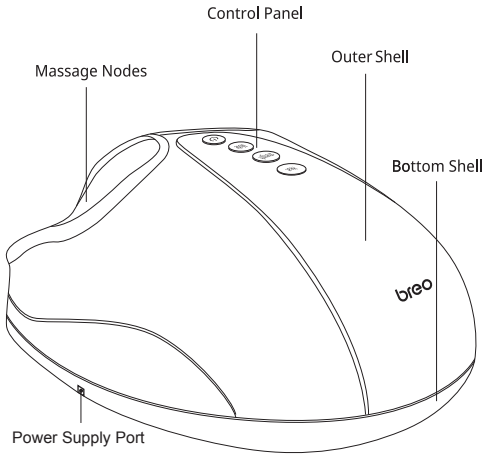Control Panel

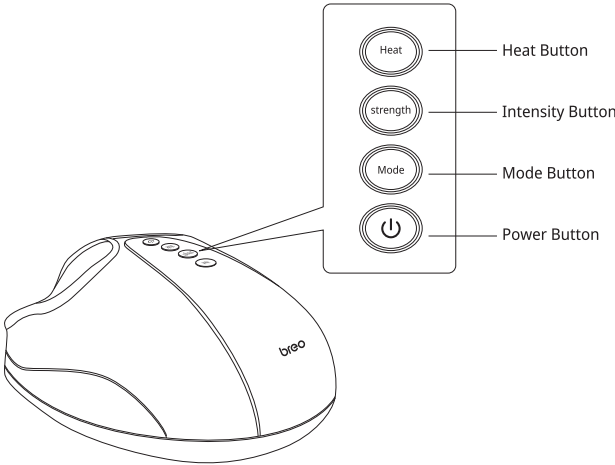Operation Guidance

- 1. Power on
  - ①Connect the power cord into the power supply port of the foot massager.
  - ② Plug the charger into the power socket.
  - ③Short press the power button to power on /off the foot massager.
  - ④The appliance is only to be used with the power adaptor provided with the appliance.
- 2. Four Buttons:  
Power On/Off Button, Mode Button, Intensity Button, Heat Button.
- 3. Power On/Off Button: When the device is turned off, short press the button to start up. Mode 1 and Intensity 1 as default.
  - ①The heating automatically turns on; the Power Button flashes red; Mode Button, Intensity Button, and Heat Button flash white.
  - ②Short press the power button to turn off the device.

- 4. Mode Button: Short press the Mode Button to switch among three massage modes.  
Mode 1 (White Light): low-speed rotation with interval air pressure massage.  
Mode 2 (Green Light): medium-speed rotation with interval air pressure massage.  
Mode 3 (Red Light): fast-speed rotation with interval air pressure massage.
- 5. Intensity Button: Short press the Intensity Button to cycle three different intensity levels.  
Intensity 1 (flashes white) -> Intensity 2 (flashes green) -> Intensity 3 (flashed red) -> Intensity 1 (flashes white)
- 6. Heat Button: The device is turned on with heat function on (flashes white), short press the Heat Button to deactivate the function.
- 7. Auto-off: The device will turn off automatically after 15 minutes of use.

Symbol description

	Medical device		Type of protection against electric shock: Class II
	Applied part of type BF		Refer to instruction manual.
	Batch code		Manufacturer information: The manufacturer name and address
	Disposal in accordance with Directive 2002/96/EC (WEEE)		Keep Dry
	Fragile handle with care		This way up
	Stacking limit by number, the maximum stacking is 4.		This device has not been tested for use in an MR environment and should not be used exposed to MR environments while patients are wearing the device. Keep it outside the MRI scanner room.

Cleaning and Maintenance

Cleaning

- 1.Before cleaning, turn off the power and unplug the power adapter from the socket. Use a soft cloth with alcohol to clean the surface of the product.
- 2.Do not use corrosive cleaners such as benzene or diluents to clean the product, otherwise it will cause discoloration, surface aging, structural damage and other conditions.

Storage

- 1.If you do not use it for a long time, please cover it with a storage bag, and place the product in a dry and ventilated room to avoid sun exposure.
- 2.This product contains batteries, do not approach high temperature or put in the fire.

Disposal

 At the end of the product life cycle, do not throw this product into the normal household garbage, but bring it to a collection point for the recycling of electronic equipment and battery. Waste Electrical and Electronic Equipment may have potentially harmful effects on the environment. Improper disposal may lead to the accumulation of harmful toxins in air, water and soil, which is harmful to human health. you have obligation to dispose of device correctly. Consult your municipal authority for information about disposal.

Customer Support

If you need help regarding product issues, after-sales services, or any other relevant inquiries, please don't hesitate to contact us:

Email: [support@breocare.com](mailto:support@breocare.com)

1-866-546-3587 (United States)

Facebook support: <https://www.facebook.com/breowellness>

For more information regarding our products, please visit [us.breo.com](https://us.breo.com)

Find Breo massagers on Social

Get connected to the community and stay updated with [#breocare](#) and [#simplyrelax](#)

 [@breocare: https://www.instagram.com/breocare/](https://www.instagram.com/breocare/)

 [@breomassagers: https://www.facebook.com/breowellness](https://www.facebook.com/breowellness)

This device complies with Part 15 of the FCC Rules. Operation is subject to the following two conditions: (1) this device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.

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