

SIMPLY RELAX

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breo
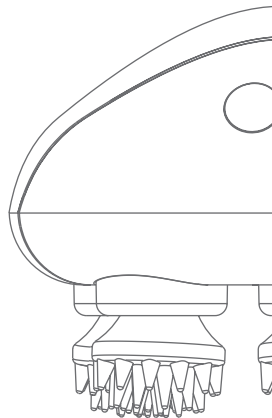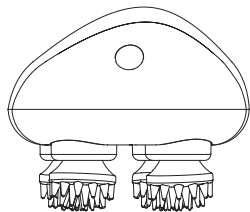
## MENU

Caution	4
Product Configuration	5
Product Characteristics	6
Product Functions	7
Troubleshooting	8
Protection	9
Product Information	10

Thank you for your purchasing Breo!

Please read this user's manual carefully prior to use and, keep it properly for future reference.


## Caution

Note: To ensure safety in operation, please read the Use Instruction Manual carefully prior to operation. The Company should not be liable for the consequences of improper use or non-compliance with the instructions for the product.

1. Do not use the product when driving a vehicle.
2. Prior to operation, remove the jewelry or articles that may affect the massage from the body part to receive massage.
3. Please keep the product away from water source, source of ignition and corrosive environment.
4. Observe safety regulations when you use the product in a specific environment.
5. Do not use the massager on body parts subjected to surgical operation.
6. Do not use the massager on body parts suffering from skin injury, swelling, inflammation or blood stasis.
7. Those with heart disease or a history of heart disease, acute disease, skin disease, malignant tumor or physical abnormality should not use the product.
8. Those with blood disease showing bleeding tendency, osteoporosis, soft tissue disease and apoplexy should not use the product.
9. Those with leaden feet, physical and mental illness, or mental retardation should not use the product, and should follow the doctor's advice if they need to use it.
10. If you have any discomfort during use, stop using the product.
11. Pregnant women and menstrual women should not use the massager.
12. Children under 6 years of age are prohibited from using the product, and children aged 6-12 can use the product only when they are subject to guidance and care from their parents or supervisors.
13. Do not scratch or stab the product using a sharp object.
14. This product is a domestic electronic health-care product. The product is not designed to work as a medical device, and it does not apply in diagnosis, treatment and other commercial uses.
15. Never fail to charge the product using charging equipment provided by the Company. Pay attention to the electrical safety.
16. Once any exception occurs during operation, power off the product immediately. Please do not disassemble, reassemble and repair the product without consent.

## Product Configuration

Massager, charging cable, storage bag, and instruction manual
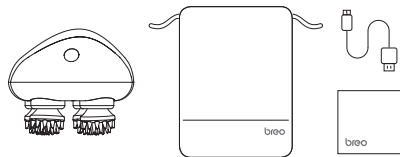
## Product Characteristics

Adopting innovative  
grasping

IPX7 waterproof

3 optional modes

Light weight & small size

## Product Functions

- 1. Start the machine:** Under the shutdown state, press the ON/OFF key, the product starts running;
- 2. Mode selection:** Short press the Mode key to switch the product between three modes.  
Mode 1: The massage head rotates fast forward, the white indicator flashes once and then keeps on, and Mode 1 is the default mode;  
Mode 2: The massage head rotates forward slowly, the white indicator flashes twice and then keeps on;  
Mode 3: The massage head rotates forward slowly and circularly, the white indicator flashes for three times and then keeps on;
- 3. Shutdown:** Long press the "ON/OFF" key to shut down the product;  
The product automatically shuts down after working for 10 minutes by default;
- 4. Electric consumption:** The indicator flashes when the product battery is low;
- 5. Charging:** The charging time is about 2~3h, during charging, the indicator flashes; when it is fully charged, the indicator keeps on;  
The product enables the self-protection program under charging state, and stops running.

## Troubleshooting

If a fault occurs during operation, do check it. In case the fault cannot be eliminated, or the product still fails to work satisfactorily, please turn off the massager and ask a professional technician to repair it. Please do not disassemble the product without consent!

Basic fault	Possible cause	Basic troubleshooting method
Failure to start up	Check whether the battery level is low	If yes, please charge in time
Low massage strength	Check whether the battery level is low	If yes, please charge in time

## Protection

### Storage

To use the machine when a storage period of over 3 months expires, never fail to perform an open-case inspection and recharge it. Never allow the machine to abut against the floor, wall, and roof when you place it. The room should be well ventilated, without intense sunlight exposure and corrosive gas erosion.

### Transportation

The product package with simple anti-vibration functions applies in air, railway, highway and ship transportation, and should be protected from rain, snow, splash, inversion and collision.

### Power Check

Check the power supply. If the supply voltage exceeds the specified range, the machine should not be turned on.

### Maintenance

Damp-proof, anti-mouse and anti-pollution measures should be taken. The machine shall be shut down before it is cleaned. Place the product in a safe place after use, and do not put any heavy object on it. The machine shall be repaired by professional technician in case any fault occurs. The non-professional are not allowed to disassemble and reassemble the product without consent.

### Abandonment and Recovery

This product should not be discarded along with household waste. Please follow local environmental laws and guidelines when you decide to dispose the product.

## Product Information

Product name: Scalp Massager

Model: Scalp 2M

Voltage: 5V==

Power: 5W

Battery: 1150mAh

Net weight: 215g


Dimensions: L98×W91×H77(mm)

It must only be supplied at safety extra low voltage corresponding to the marking on the appliance.

This appliance contains batteries that are only replaceable by skilled persons.

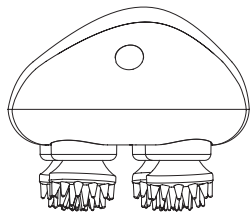This appliance contains batteries that are non-replaceable.

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