breo iPalm 520E

USER
MANUAL

breo


Thank you for purchasing Breo iPalm 520E.

Please read the manual carefully before using the product, and keep it in a safe place for future reference.

Inside the Package

Main Unit, Charging Cable, User Manual

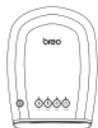Main Unit

Charging Cable


User Manual

A. Operating Buttons

From left to right: ①Air Pressure ②Heat ③Time ④On/Off

B. Air Pressure Reduction Valve

C. Power Port

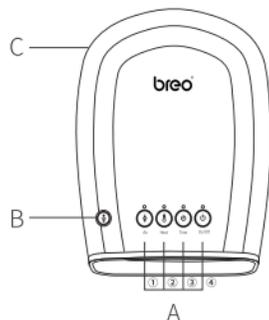01

Preparation

Operations

Before use, check the condition of the machine and accessories:

1. Battery power is sufficient.
2. No damage on adapter and charging cable.

Charging Connection

1. Connect the adapter with the charging cable (USB end).
2. Insert the other end of the charging cable into the main unit's DC port.
3. Connect the plug of the adapter to the AC power socket.

Note:

- The red LED light flashes when charging, and the light turns green when the device fully charged.
- The device can not work when charging.

Operation

Use

- Place hand (palm down) into the massager.
- Turn the power on, and set air pressure & heat temperature & time for preference.

02

Function Buttons

On/Off

Power On—Press and hold the On/Off Button until the red LED light turns on.

Power Off—Press and hold the On/Off Button until the red LED light turns off.

Note: After turning off the unit, allow approximately 30 seconds to fully release air pressure.

Time

Press the Time Button to switch massage duration (5min/10min/15min). The LED light for Time flashes and the massage duration circulates among these three modes:

1. 5 minutes (default): the LED light flashes once.
2. 10 minutes: the LED light flashes twice.
3. 15 minutes: the LED light flashes three times.

Heat

Press the Heat Button to select temperature (High-Off-Low).

The LED light for Heat flashes and the temperature circulates among these three modes:

1. High (default): the LED light flashes twice
2. Off: the LED light turns off.
3. Low: the LED light turns on.

Air

Press the Air Button to select pressure level (Default-Standard-Powerful-Soft)

The LED light for Air Pressure flashes and the setting circulates among these four modes:

1. Circulating from low to high (default): the LED light flashes once.
2. Standard: the LED light flashes twice.
3. Powerful: the LED light flashes three times.
4. Soft: the LED light flashes four times.

Air Pressure Reduction Valve Button

If you feel uncomfortable during the massage session, press the Pressure Reduction Valve (red button) to release air, take your hand out and turn off the device.

Preset Mode

The device will operate automatically when it is turned on. The preset settings are circulating air pressure strength (from soft to powerful), high-temperature heat compression, 5 minutes.

Customized Mode

Press related function buttons to select air pressure level, hot compress level and massage time.

Precautions

To ensure safe and correct use of the product, please read the Operation Manual in detail before use. Our company will assume no liability for any consequences arising from use of the product improperly or without observing the instructions.

1. Do not use the product during driving.
2. Prior to use, please remove any jewelry or object that may affect massage from the massage position.
3. Do not use the product in bathroom or other damp environments.
4. Keep the product away from water and fire sources and corrosive environments.
5. Please observe safety regulations when using the product in specific environments.
6. Do not use the product if the massage position has undergone operation.
7. Do not use the product if the massage position has skin injury, redness and inflammation, congestion or other discomfort symptoms.
8. Persons with heart disease or history of heart disease, or with acute disease, skin disease, malignant tumor or physical abnormality should not use the product.
9. Persons with bleeding tendency or other blood diseases, or with osteoporosis, soft-tissue disorders or apoplexia should not use the product.
10. Persons with slow action, psychosomatic disorders or disturbance of

intelligence should not use the product; please follow the doctor's advice if any of these persons intends to use the product.

11. Stop using the product in case of any discomfort during use.
12. Women who are pregnant or in menstrual period should not use the massager.
13. Children under 8 years of age are forbidden to use the product; children between 8-12 years of age must use the product under the guidance and care of their parents or other persons with supervision capability.
14. Use alcohol cotton to clean the external appearance; avoid using water for cleaning.
15. Do not use any sharp objects to scratch or stab the product.
16. The appliance has a heated surface. Persons insensitive to heat must be careful when using the appliance.

Battery removal steps

Operation steps for battery removal

1. Cut off the external power supply;
2. Disconnect the power supply to the battery;
3. Remove the battery;
4. Place the waste battery at a fixed recycling station;
5. Please observe local environmental protection laws and waste disposal guidelines to protect the environment.

Notes

- Remove the battery out of the appliance before scrapping the appliance;
- The appliance must be disconnected from power supply before removing the battery;
- The battery should be disposed of safely.

Except professionals, no one shall attempt to remove the battery!

Troubleshooting

If a fault occurs during the operation of the device, please check the following precaution.

If the error cannot be fixed or the device still does not work properly, please cut off the power and contact Breo to repair the device. Do not disassemble or try to repair the device yourself.

Problem	Possible Reason	Possible Solution
Power-On failure	Batteries are low / out	Charge the battery
Weak heat or air pressure	Batteries are low / out	Charge the battery

If you have any issues, questions or feedback, please email us at support@breo-care.com or call us at 866-546-3587. We are always happy to help.

Protection

Storage

When the storage period exceeds 3 months, unpack the product, and check and recharge it before use. Do not place the product tightly close to floor, wall and ceiling.

The indoor environment should be well-ventilated; keep away from strong sunlight and corrosive gas.

Transport

The product package is provided with simple shockproof measures, therefore is suitable for transport by air, rail, road and sea. Keep away from rain and snow, and avoid inversion and collision.

Check of power supply

Check the power supply; do not power on the product if the supply voltage is beyond the specified range.

Maintenance

Take measures to resist moisture, rats and pollution. Cleaning should be carried out when the product is powered off. After use, put the device in a safe place; do not place any heavy objects on the device. In case of any failure, please ask a professional technician to repair it. Do not attempt to disassemble the product.

Disposal and Recycle

The product should not be abandoned along with domestic wastes. Dispose of the product according to local environmental laws and guidelines.

Precautions

1. Please remember to respect the local regulations: hand in the non-working electrical equipments to an appropriate waste disposal center. The packaging material is recyclable. Dispose of the packaging in an environmentally friendly manner and make it available for the recyclable material collection-service.
2. This appliance can be used by children aged from 8 years and above and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge if they have been given supervision or instruction concerning use of the appliance in a safe way and understand the hazards involved.
3. Children shall not play with the appliance. Cleaning and user maintenance shall not be made by children without supervision.

Product information

Name	Hand Massager
Model	iPalm 520E
voltage	5V ---
power	5W
Battery	2300mAh
weight	1060g
Dimension	L265×W201×H90(mm)
Ideal environment	Temperature: 5°C-40°C Relative Humidity: ≤80% Atmospheric pressure: 860hpa~1060hpa

Made in China
Shenzhen Breo Technology Co., LTD
19/F, Hisense Southern Building, No.1777 Chuangye Road, Nanshan
District, Shenzhen 518054, China

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