

---

Thank you for purchasing Breo iNeck2  
Please read this manual carefully before using  
and keep it properly for future reference.

---

## Contents

---


Inside the Package	3
Function Deployment	4
Preparation	5
Operation	6
Precautions	7
Troubleshooting & Technical Parameters	8
Storage & Transportation	9



## Inside the Package

---

Your massager comes with a Master Unit, Power Charger, User Manual, Portable Bag.

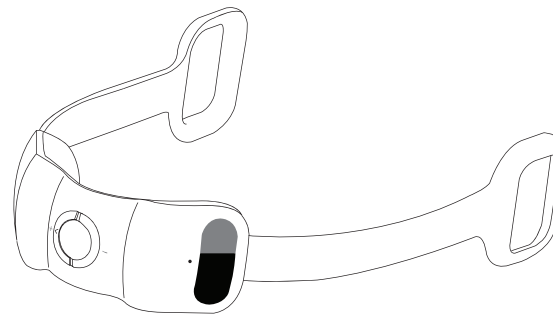

Master Unit




Power Charger



User Manual



## Function Deployment

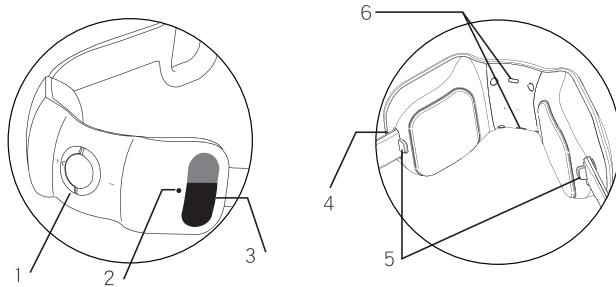

- [ 1 ]** drawstring width knob
- [ 2 ]** LED indicator
- [ 3 ]** function Operating Buttons (including LCD display, mode selection key 、 On/Off key)
- [ 4 ]** Power charging jack
- [ 5 ]** Remove the chain with key
- [ 6 ]** two infrared lamps

## Preparation

### Preparation

1. Insert the plug of power charger into power charging jack.
2. The indicator of power on the LCD screen flashes during charging,the flash stops when fully charged.
3. In the process of charging, the product does not support operating (if the product is charged while operating, it will start the self-protection program and stop operating automatically)
4. Charging time is 5-6 hours.

### Wear

Pull belt into belt hole, please use nails on massage head width adjustment knobs concave, button out V half arc, put on the main unit to adjust drawstring width knob till a desired position got,please aim the main unit right at neck, hold drawstring with two hands,

### Power-on

If to short press Power On/Off key, the product starts working, and LED indicator lights in red. The default time is 15 minute, and the default massage mode is mode 1. To change massage mode, please operate mode repeat key, and the corresponding massage mode shall be displayed on LED display.

### Manual mode

Mode selection: four massage modes: if to short press mode selection key, massage mode will be change in turn.

**Mode1** one group of massage head runs forward and one group of massage head runs reversely;

**Mode2** one group of massage head runs reversely and one group of massage head runs forward;

**Mode3** two groups of massage head run forward (with infrared magnetotherapy function);


**Mode4** two groups of massage head run reversely (with infrared magnetotherapy function) ;

## Operation

### Massage strength adjustment

Intellectual adjustment of massage strength: when the massage head touches articles, massage strength will automatically increase; when the massage head fails to touch articles, massage speed will automatically decrease.

- [1]** symbol of battery level
- [2]** Massage mode including mode 1, 2, 3 and 4 displayed
- [3]** Massage time including massage time 5 minutes, 10 minutes and 15 minutes displayed
- [4]** Symbol of repeat
- [5]** Symbol of time
- [6]** Mode select key
- [7]** On/Off Key

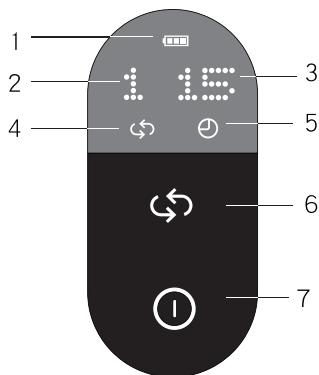

### Power-off

If to long press Power On/Off key, the products shall stop running, and LED turns off. (When fixed time is up, the product will automatically turn off).

If to short press Power On/Off key, running time will switch to 5min, 10min and 15min, which will show on LCD screen.

Please do read this User's Manual carefully for correct and safe use of this massager.

1. Do not use the neck massager while driving.
2. Do not use the massager at the place with excessively high temperature and moisture.
3. People suffering from swelling, inflammation or blood stasis in the neck, or have had neck surgery shouldn't use the massager.
4. Patients with apoplexy, high blood pressure and hypoglycemia are prohibited to use the product.
5. Those slow in action, discomfort in body and mind or the mentally handicapped are prohibited using this product.
6. Please avoid using this neck massager after drinking and when body unwell.
7. This massager is designed for household health care, and does not apply to the diagnosis, treatment and other commercial uses.
8. Please wipe the massager with sterile cotton moistened with neutral detergent. Do not clean it with water. Switch off the device before cleaning.
9. The massager has hot compress function, users insensitive with temperature sensing should be careful in use.
10. Noise during massage is normal.
11. The massager is intellectually controlled. In case of over large of strength, it starts machine protection mechanism, and it is normal to stop running.
12. Check the power supply and donot switch on the device if the voltage exceeds the rated range of the device.
13. This appliance is not intended for use by persons (including children) with reduced physical, sensory or mental capabilities, or lack of experience and knowledge, unless

## Precautions

---

they have been given supervision or instruction concerning use of the appliance by a person responsible for their safety.

14. it must only be supplied at safety extra low voltage corresponding to the marking on the appliance.
15. This appliance contains batteries that are only replaceable by skilled persons.
16. This appliance contains batteries that are non-replaceable.

## Troubleshooting & Technical Parameters

---

If there is any problem when you use the device, please follow the instructions below to check. If the problem still exists or the device can't work normally, please turn off the power supply and contact local service center. Do not disassemble or repair the device yourself.

Fault	Possible Reasons	Solutions
Power-on fail	Whether the battery is lows	If yes, Please charge at once
The vibration is over low.	Whether the battery is low	If yes, Please charge at once

## Troubleshooting & Technical Parameters

---

### Technical parameters

Power	5V == 1A
Ideal environment	Magnetic field strength:50mT ~ 250mT Fomentation:30 ℃ ~ 45 ℃
Dimension	Master unit : Length 220mm×Width 105mm× Height 120mm
N.W.	855g

## Storage & Transportation

---

### Storage

In case the storage period exceeds 3 months, please check it and power it on for half an hour. Do not store close to floor, wall or roof.Do not put heary objects on the device.Keep good ventilation condition and away from direct sunshine and corrosive gases.

### Transportation

The product has simple shockproof design and can be transported by air, railway, highway or ship; avoid rain, snow, impact, and always keep the right side up.