breo.iPalm 520 Pro

USER
MANUAL


breo

Thank you for purchasing Breo iPalm 520 Pro.

Please read the manual carefully before using the product, and keep it in a safe place for future reference.

Inside the Package

Main Unit, Charging Cable, User Manual

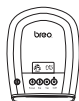Main Unit

Charging Cable

User Manual


A. Operating Buttons

From left to right: ① Air Pressure ② Heat ③ Time ④ On/Off

B. Air Pressure Reduction Valve

C. LCD Screen

D. Power Port

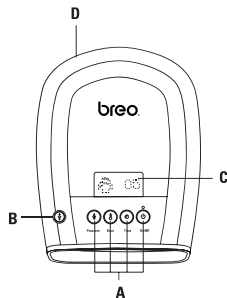01

Preparation

Operations

Prior to use, check the condition of the machine and accessories:

- Battery power is sufficient
- No damage on adapter and charging cable

Charging Connection

- Connect the adapter with the charging cable (USB end).
 - Insert the other end of charging cable into the main unit's DC port.
 - Connect the plug of adapter to AC power socket.
- Please check if the device and accessories are complete.
Please check if the power adapter is damaged.

Operation

Use

- Place hand (palm down) into the massager.
- Turn the power on, select massage mode, and set air pressure & heat at your preference.

Manual Mode

Power On/Off

Power On—Press and hold On/Off Button until the light turns on. Preset mode is air pressure (circulating from low to high) / with heat / 5 minutes duration.

02

Press related function buttons to switch between air pressure and hot compress modes.

Power Off—Press and hold ON/Off Button

*Note: After turning off the unit, allow approximately 40 seconds for air pressure to be released. During this period, the unit will not start when the On/Off button is pressed.

Time—Press the Time Button to switch massage duration (5min/10min/15min)

Power On/off

Power on—Press Button

Preset default is air pressure (Auto: Cycles through Low-Med-High).
high heat, 5 minutes duration

Manually operate related function buttons to change setting for air pressure heat compression and time.

Power Off—Press Button

Note: It will take a few seconds before the LCD display turns.

Manual mode

Air Pressure —Press Button to change massage modes.

Auto (default)—Cycles through Low-Medium-High

*Lines on the sides of hand on LCD display are flashing. These lines change from 1 line on each side for Low, to 2 lines for Medium, then 3 lines for High. The cycle starts over from low to high for the entire massage duration or until the Pressure button is pressed to manually change the massage mode.

Pressure 1—Low

*Lines on the sides of hand on LCD display are fixed—1 line on each side.

Pressure 2—Medium

*Lines on the sides of hand on LCD display are fixed—2 line on each side.

Pressure 3—High

*Lines on the sides of hand on LCD display are fixed—3 line on each side.

*Important Tip: When you feel the air pressure is too much and/or uncomfortable, press the Pressure Reduction Valve (red button) to release air.

Precautions

To ensure safe and correct use of the product, please read the Operation Manual in detail before use. Our company will assume no liability for any consequences arising from use of the product improperly or without observing the instructions.

1. Do not use the product during driving.
2. Prior to use, please remove any jewelry or object that may affect massage from the massage position.
3. Do not use the product in bathroom or other damp environments.
4. Keep the product away from water and fire sources and corrosive environments.
5. Please observe safety regulations when using the product in specific environments.
6. Do not use the product if the massage position has undergone operation.
7. Do not use the product if the massage position has skin injury, redness and inflammation, congestion or other discomfort symptoms.
8. Persons with heart disease or history of heart disease, or with acute disease, skin disease, malignant tumor or physical abnormality should not use the product.
9. Persons with bleeding tendency or other blood diseases, or with osteoporosis, soft-tissue disorders or apoplexia should not use the product.
10. Persons with slow action, psychosomatic disorders or disturbance of

intelligence should not use the product; please follow the doctor's advice if any of these persons intends to use the product.

11. Stop using the product in case of any discomfort during use.
12. Women who are pregnant or in menstrual period should not use the massager.
13. Children under 8 years of age are forbidden to use the product; children between 8-12 years of age must use the product under the guidance and care of their parents or other persons with supervision capability.
14. Use alcohol cotton to clean the external appearance; avoid using water for cleaning.
15. Do not use any sharp objects to scratch or stab the product..
16. The appliance has a heated surface. Persons insensitive to heat must be careful when using the appliance.

Battery removal steps

Operation steps for battery removal

1. Cut off the external power supply;
2. Disconnect the power supply to the battery;
3. Remove the battery;
4. Place the waste battery at a fixed recycling station;
5. Please observe local environmental protection laws and waste disposal guidelines to protect the environment.

Notes

- Remove the battery out of the appliance before scrapping the appliance;
- The appliance must be disconnected from power supply before removing the battery;
- The battery should be disposed of safely.

Except professionals, no one shall attempt to remove the battery!

Troubleshooting

If a fault occurs during the operation of the device, please check the following precaution.

If the error cannot be fixed or the device still does not work properly, please cut off the power and contact Breo to repair the device. Do not disassemble or try to repair the device yourself.

Problem	Possible Reason	Possible Solution
Power-On failure	Batteries are low / out	Charge the battery
Weak heat or air pressure	Batteries are low / out	Charge the battery

If you have any issues, questions or feedback, please email us at support@breo-care.com or call us at 866-546-3587. We are always happy to help.

Protection

Storage

When the storage period exceeds 3 months, unpack the product, and check and recharge it before use. Do not place the product tightly close to floor, wall and ceiling.

The indoor environment should be well-ventilated; keep away from strong sunlight and corrosive gas.

Transport

The product package is provided with simple shockproof measures, therefore is suitable for transport by air, rail, road and sea. Keep away from rain and snow, and avoid inversion and collision.

Check of power supply

Check the power supply; do not power on the product if the supply voltage is beyond the specified range.

Maintenance

Take measures to resist moisture, rats and pollution. Cleaning should be carried out when the product is powered off. After use, put the device in a safe place; do not place any heavy objects on the device. In case of any failure, please ask a professional technician to repair it. Do not attempt to disassemble the product.

Disposal and Recycle

The product should not be abandoned along with domestic wastes. Dispose of the product according to local environmental laws and guidelines.

Precautions

1. Please remember to respect the local regulations: hand in the non-working electrical equipments to an appropriate waste disposal center. The packaging material is recyclable. Dispose of the packaging in an environmentally friendly manner and make it available for the recyclable material collection-service.
2. This appliance can be used by children aged from 8 years and above and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge if they have been given supervision or instruction concerning use of the appliance in a safe way and understand the hazards involved.
- 3.Children shall not play with the appliance. Cleaning and user maintenance shall not be made by children without supervision.

Product information

Name	Hand Massager
Model	iPalm 520 Pro
voltage	5V---
power	5W
Battery	2300mAh
weight	1060g
Dimension	L265×W201×H90(mm)
Ideal environment	Temperature: 5°C-40°C Relative Humidity: ≤80% Atmospheric pressure: 860hpa-1060hpa

Made in China
Shenzhen Breo Technology Co., LTD
19/F, Hisense Southern Building, No.1777 Chuangye Road, Nanshan District, Shenzhen 518054, China

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