


NURSAL®

PORTABLE ORAL IRRIGATOR

Model: DC5121

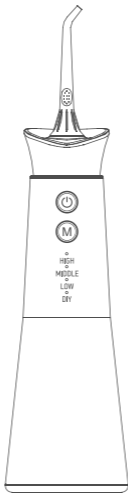

USER MANUAL

CONTENTS

Important safety measures	01
Product description	05
Product accessories	06
Product specifications	07
Using the appliance	08
Water pressure adjustment	10
Tips	12
Unit disassembly	13
Battery maintenance	14
Troubleshooting	15
After-sales service	16

IMPORTANT SAFETY MEASURES

Basic safety precautions should be followed when using electrical products, especially when there are children in the house.

Please read all the outlined instructions before the first use.

IMPORTANT!

To reduce the risk of electric shock:

DO NOT handle the charger with wet hands.

DO NOT immerse the charger in water or any other liquids.

DO NOT place or store the product in a place where it may fall into the bathtub, shower, or sink.

When charging, if the product falls into the water, please unplug the unit immediately before picking it up.

Dry the product before recharging or using it again.

Please check whether the charger cable is damaged before the first use.

WARNING:

In order to reduce the risk of electric shock, fire, injury, or damage:

Do not attempt to open or disassemble your DC5121 Portable Oral Irrigator, it will void the warranty.

Do not plug the unit into an adapter that does not meet the specified voltage. If the product power cord or plug is damaged, or the product was dropped, damaged, or falls into water, do not use it.

Please contact your local dealer.

This appliance is designed for use only in the oral cavity. DO NOT direct the water flow under the tongue, into the eyes, ears, nose, or other minor organs. This product is capable of producing pressure that may cause serious damage to these organs. Please refer to the operating instructions for proper use.

Magnetic fields: Pacemakers and similar implanted medical devices may be affected by magnetic fields. A small amount of bleeding is normal when using the unit for the first time.

If the gums bleeding is severe, stop using the oral irrigator and contact your dentist.

Fill the water tank with lukewarm water or other solution recommended by your dentist.

NEVER fill the water tank with hot water! Sensitive gums can be easily irritated.

Do not insert nor drop any foreign objects into the inside through the openings.

Keep the charger and its wires away from heated surfaces and heat-generating appliances. Do not use outdoors or operate where aerosol (spray) products are being used or where oxygen is being administered.

Do not use iodine, salt water (saltwater solution) or waterinsoluble concentrated essential oils in this product. The use of such substances will reduce unit performance and shorten its lifespan. If these items are used, it may invalidate the warranty.

Do not use while wearing oral accessories, and please remove them before use.

Do not use an oral irrigator if there is an open wound on your tongue or mouth.

If you are suffering from tooth or gum disease or receiving oral treatment, consult a dentist before using this oral irrigator or any other oral hygiene auxiliary chemical.

When children (over eight years old) and special groups of people do use this product, they must be guided and supervised closely.

Make sure this product is used properly.


Prevent children from playing with this electrical product.

There are no user-maintainable parts in this appliance, and no maintenance is required.


Only fill the water tank with lukewarm water (below 30°C) or other solutions recommended by your dentist.

Please use this product in accordance with the user guide or recommendations from your dentist.


PRODUCT DESCRIPTION

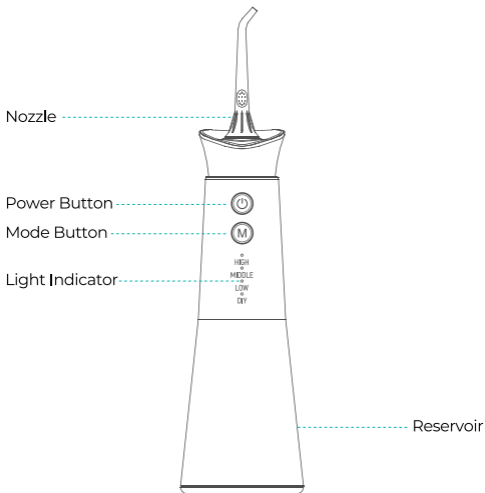

PRODUCT ACCESSORIES

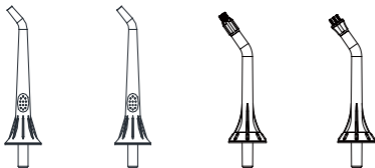

Nozzle x4

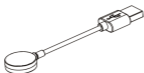

Charging cable x1

Replaceable nozzles/accessories can be purchased from the vendor.

PRODUCT SPECIFICATIONS

Product Name		NURSAL Protable
Oral irrigator	Product number	DC5121
	Water tank capacity	300ML
	Rated working voltage	3.7V ===
	Rated power	5W
	Waterproof strength	IPX8
	Charging time	Requires more >4 hours to fully charge
	Input parameters	5V 1A
	Origin	China
	Executive standard	GB4706.1-2005 GB4706.59-2008

USING THE APPLIANCE


CHARGING THE UNIT

VERY IMPORTANT!


Please fully charge the unit before using it for the first time.

Plug the USB cable into the USB port on the back of your portable oral irrigator and turn it on.

The product cannot be used while charging.

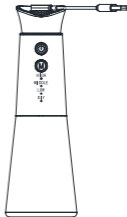

Open the reservoir cap to fill it with the lukewarm water
Note: Please make sure that the device is turned off when filling the reservoir.

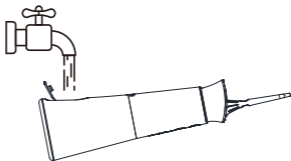

ATTACH/DETACH THE NOZZLE

Nozzle loading and unloading


The nozzle is assembled by inserting the large hole of the nozzle into the water outlet hole at the top of the body (the nozzle is magnetically fastened to the body).

To remove the nozzle, pull the nozzle upward by hand to remove it.

After the nozzle is installed, the nozzle angle can be rotated.

After the nozzle is installed, the nozzle must not come loose when using the highest water pressure.

The installed nozzle must be a DC5121 replacement nozzle.

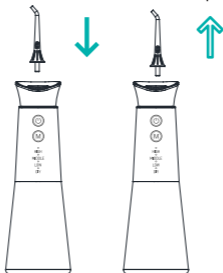

WATER PRESSURE ADJUSTMENT

In the off state, press the "M" button continuously, the mode in LOW- MEDIUM- HIGH-DIY (Clean)-DIY (Soft)-DIY (Pulse)-DIY (RINSE) to switch between, select the appropriate mode and then turn on, the first time you use this product, please start from LOW gear, after the oral adaptation can gradually adjust the high level.

As each person's gum tolerance is different, there may be gum bleeding, this is a normal phenomenon, please feel free to use! If there is a lot of bleeding, please stop using and consult your dentist.

Long press the power button to turn on/off, any mode continuous working time is 90S, This product has memory function, when you use it next time, it will continue the mode you have used last time.



DIY(Clean)

DIY(Soft)

DIY(Pulse)

DIY(RINSE)

TIPS

Lean over the sink, press the mode button to adjust the water pressure to the desired level, place the nozzle head into the mouth. Hold the unit in an upright position with the nozzle tip directed toward your teeth.

Press and hold the power button for 2 seconds to turn the unit on.



Align the nozzle head with the gums at an angle of 90 degrees so that the bristles are gently touching the teeth. Start flushing from the molar teeth (back teeth) towards the front teeth. Gently move the nozzle along the gums with a short pause between the teeth to allow the water to pass between the teeth. Repeat the same procedure on each side of the teeth.




UNIT DISASSEMBLE

Turn off the unit.


Open the reservoir cap with the hole facing down as shown in the figure below.

Pour the liquid out, turn the unit on to drain the residual liquid until no more water is left in the reservoir.

Turn off the device.

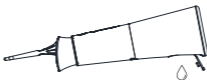

Grab the unit and rotate the upper unit body part counterclockwise to remove the reservoir as shown in the figures below.

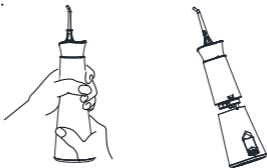

Note:

There may be moisture inside the unit before the first use (in the reservoir and suction pipe), which is the residue of the distilled water used for product testing.

BATTERY MAINTENANCE

Battery maintenance

The device must be fully charged before the first use, charging takes about 4 hours;

Plug the USB cable into the charging socket and connect it to a power source;

The battery needs to be fully charged before storing the unit for a long time period with no use;

In order to extend the battery lifespan, please charge the battery immediately when it is low.

TROUBLESHOOTING

Problem	Cause	Solution
Leakage between the nozzle and the nozzle insertion hole	The nozzle does not fit well with the handle or the seal ring is damaged	Insert the nozzle firmly into the nozzle insertion hole until it clicks into place; Contact the vendor in order to repair your unit
Water is leaking through the reservoir cap	The seal ring is damaged, or there are impurities on the inlet side	Contact the vendor for repair
Water tank is leaking	The seal ring is damaged	Replace the seal ring, contact the vendor for repair
Low water pressure	The battery is low; An adapter that does not meet the specified voltage is in use	Please fully charge the unit
Unit does not work	The battery is low	Please fully charge the unit; Contact the vendor for repair
No water flow	First use or not used for a long time	Fill the reservoir with the proper amount of water, reinstall the nozzle, power on the unit

If you still can't solve the problem, please consult the vendor or go to the maintenance place designated by the NURSAL.

AFTER-SALES SERVICE

Thank you for purchasing the NURSAL Oral Irrigator DC5121.

Please pay attention to the following statements:

1. Ask the vendor to issue a valid sales certificate as the basis for the warranty.

2. According to relevant national regulations, this product implements free maintenance for product failures due to non-human factors within one year (from the date of purchase on the invoice). The warranty does not include NURSAL. If the warranty period is exceeded, it can be sent to the company's designated after-sales maintenance place for paid maintenance.

3. During the warranty period, if the following situations occur, costs and service fees will be charged as appropriate:
Consumer-caused damage due to improper use, maintenance, and storage;
Damage is caused by a repairer who does not undertake the guarantees of return, repair, and replacement;

Lack of a relevant sales voucher and the production date of the appliance has passed one year;

Damage is caused by commercial use;

Damage is caused by force majeure factors (such as lightning, earthquake, fire, etc)

4. Please fill in the following form carefully, and send it to the after-sales repair place designated by the company for repair along with the appliance that needs repairing.

WARRANTY REPAIR RECORD

Customer name		Contact number	
Product number		Product color	
Purchase date		Purchase platform	
Contact address			
Repair date			

CUSTOMER SERVICE

✉ support@nursalshop.com

☎ 1-855-666-8877



Made in China