 **Sportneer**

*Handheld Deep Tissue Massager*

*User Manual*

# 1 Description

**This device has the following benefits:**

## **1. Warm you up before exercise and recover after exercise**

Equipped with a high-speed motor and 11mm large-amplitude head attachment, this massager activates and repairs muscles with ease.

## **2. Great for you, undistruptive to others**

With only 40dB of sound output, the percussions let you relax fully without making excessive noise.

## **3. Lightweight and pocket-sized**

Tough yet lightweight. Chrome body weighs just 1.05 kg, including the battery!

## **4. User-friendly, one-button control**

Operate the one-button control with just one hand while holding the comfortable and ergonomic handle for smooth, convenient adjustment of vibration intensity.



## **5. Pain relief from all types of activities**

If you're experiencing soreness from sitting or lifting heavy objects, the Sportneer deep tissue massager can reach those deep tissue areas to relieve tension in your neck, back, shoulders and legs.

## 2 Warnings

1. Please use the provided charger to charge the device. It takes about 3~5 hours to charge the massager fully (once fully charged, the indicator light on the power adapter will change from red to green.).
2. Operating temperature is between 41-104°F.
3. Do not operate the device while it is charging.
4. Do not insert or touch the gap between the massage head and its shell with fingers or other objects to avoid damaging it.
5. Do not try to modify or disassemble the device (except swapping the massage heads) by yourself.
6. Because of the intense nature of the massager we recommend operating the device for no more than 20 minutes at a time to avoid any physical discomfort.
7. Supervise closely if used by children
8. If you feel pain or any discomfort, be sure to stop using it immediately.
9. Consult your physician prior to use if any of the following conditions apply to you:
  - Osteoporosis
  - Heart problems
  - You have an electronic medical device
  - Pregnant women
  - You have experienced physical trauma
  - You have other medical complications


### 3 Packing List

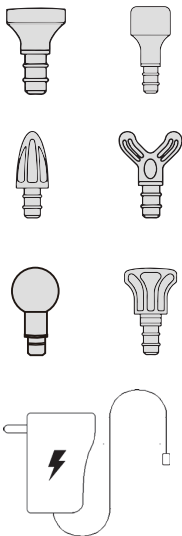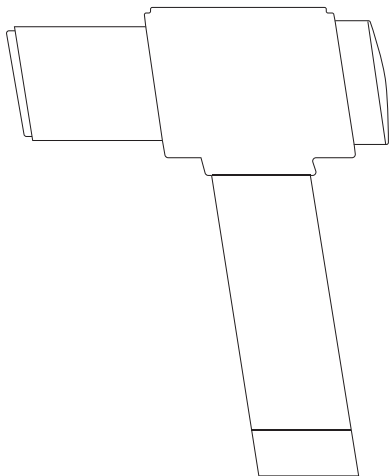
#### Technical Parameter

Model	9°
Power	12-220W
Frequency	22Hz/28Hz/33Hz/44Hz/48Hz/55Hz
Battery Capacity	3400mAh (Panasonic battery cells)/56Wh
Battery Life	1.5-6.6h
Torque Pressure	0-8.5kg/cm
Amplitude	11mm
Intelligent Mode	10min timer
Operating Temperature	41-104°F

## Product Illustration
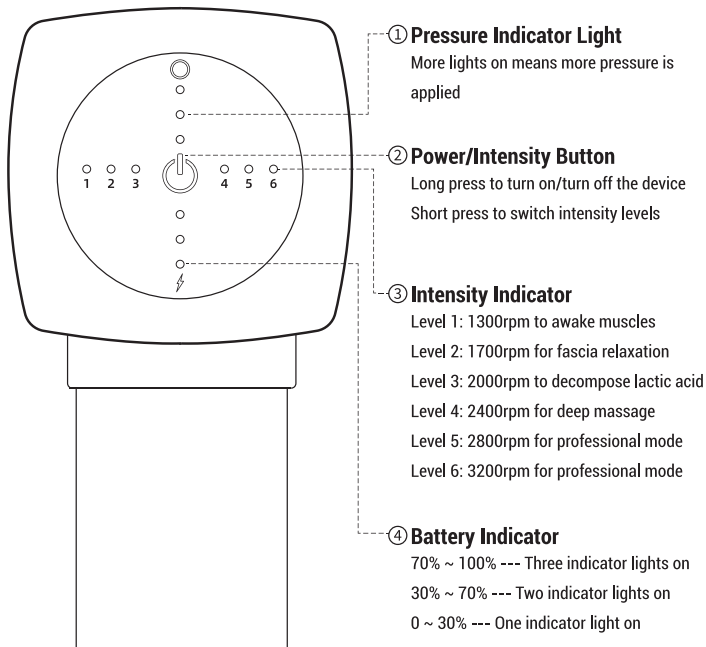
## 5 Massage Heads

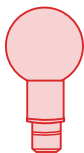**Ball Head**

For large muscle groups (quads, glutes, etc.) and extreme muscle and joint soreness

**Fork Head**

Mainly used for spine, neck, calves, and shoulder for sensitive muscles or bony areas. Intended for intensity levels 1 & 2.

**Cone Head**

Mainly used for joints, deep tissue and small muscle areas (feet, wrist, etc. ). Intended for intensity level 1 & 2 .

**Flat Head**

Ideal for general use, suitable for all body muscles

**Chrome Flat Head**

Amplified intensity, Oil applicable, Easy cleaning

**Chrome Ball Head**

Amplified intensity, Oil applicable, Easy cleaning

## 6 Charging & Operating

1. Please use the provided charger to charge the device. It takes about 3~5 hours to fully charge the massager (When fully charged, the indicator light of massager turns off, and the indicator on the power adapter will change from red to green.).
2. Pick the appropriate head attachment for your needs, then plug it into the nozzle until the two parts are fully engaged.
3. Keep the power button pressed for a moment to turn on the massager, then short press it again to switch between the different intensity levels, from level 1, 2, 3, 4, 5, 6 and then off.

## 7 Maintenance, Cleaning & Storage

This device has the 10-minute-off timer setting in case the motor inside get overheated and toasted due to overload. Long period of forcing great pressure to enhance percussion performance will often lead to overheat. If user just applies little pressure in assisting massage, you can restart this device immediately after 10-minute off.

While installing head attachments, you can twist one into the nozzle of the percussion gun instead of pushing it in. This way may prevent the gaskets on head attachments from too much fraying. We also prepared spare gaskets for you to replace when in need.

Use a slightly damp towel to clean the surface of device, then wipe it dry with a soft cloth.

Use the provided carrying case for storage or transportation. Be careful to avoid exposure to water, dust or other foreign matter that may damage it.

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