Model Item : FG850LB


Smart Body Composition Scale

User Manual

Read this manual thoroughly before using
and save for future reference

contents

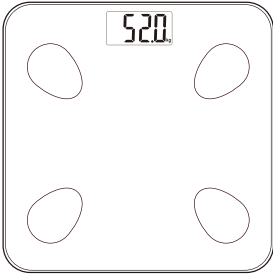1 x Scale

2 x AAA Batteries

Specifications:



Max.:180kg/400lb/28st

Min.: 6kg/13.2lb/0.93st

Division:0.1kg/0.2lb

Battery:2x1.5V AAA included

At a Glance

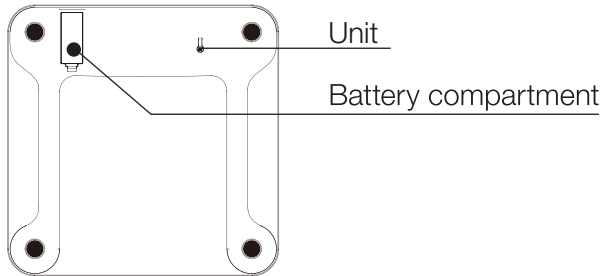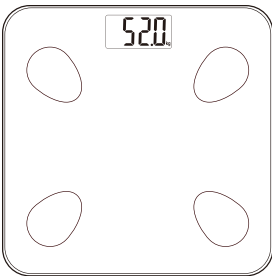Preparation

Scale Setup

- 1) Insert 2xAAA batteries into the batteries compartment.
- 2) Place your scale on a hard,level surface.

App Setup

- 1) Download the App




1. Go to IOS App store or Google play,search the keyword "Fitdays",or scan the following QR Code with your phone to find the app.
2. Download and install the "Fitdays" to your phone.


- 2) Enable Bluetooth on your smart phone.
- 3) Open the App on your phone,register your account using either your email address or your social media account.
- 4) Add your user profile by entering your Gender,Name,Birthday,Height and Target Weight.

Screen Shots

1.Register

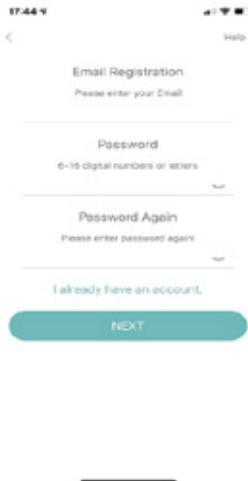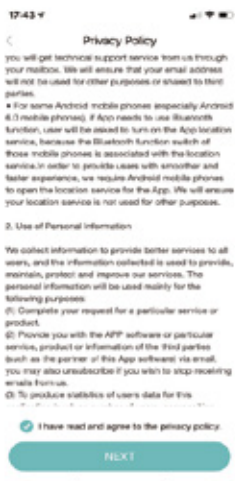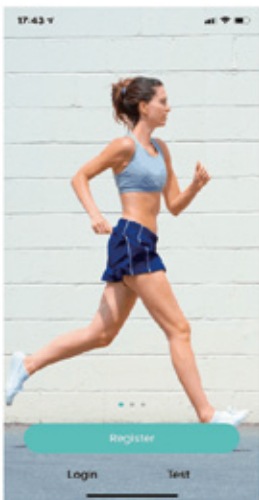1.Enter user profile

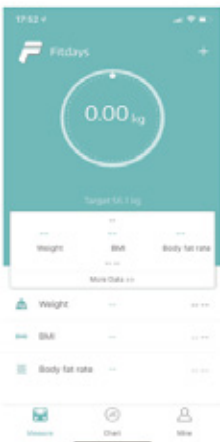Clean and Bare feet

With socks or shoes

Using the Scale

Normal Body Weighing (Without Body Composition)

- 1) Place your scale on a hard,level surface.
- 2) Step on the weighing platform with both feet and evenly distribute your weight.
- 3) Once a steady reading is taken,the weight will lock and flash to indicate the final weight reading.

Note:The scale has an auto-off feature that will automatically shut off in approximately 10 seconds without operation.

Body Composition Analyzing(With APP)


- 1) Place your scale on a hard,level surface.
- 2) Make sure the app has been properly setup and the bluetooth option is enabled on your phone.
- 3) Open the app and leave the App on at its home screen,It will automatically connect to the scale.
- 4) Step on the scale with bare feet, make sure both feet are in contact with all 4 bio-electronical impedance sensors(the metal plates).
- 5) Stand still on the scale for approx.8 seconds or until finished analyzing,
- 6) The analyzed date will be displayed and stored in your phone app.

Cautions/Warnings

- To avoid tipping, please be careful so as to not step on the edge of the scale.
- Do not hop on scales.
- Protect scale from blunt force, temperature fluctuations, and potential heat sources (e.g. stoves, heating units)
- Do not drop the scales or other objects onto it as it may damage the internal sensors.
- Do not use the scale when your feet or body is wet.
- Scale has a max capacity of 400lbs, please do not place anything on it that weights more than 400lbs.
- Place the scale on level flooring-otherwise, it may fall over.
- Clean the scale with a slightly damp cloth. DO NOT immerse the scale in water or use chemical /abrasive cleaning agents.

Battery Handling

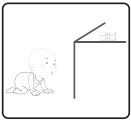-If the scale will be stored for a long period, it is recommended to remove the batteries to avoid leakage.

-Keep the batteries out of reach of children.


-Do NOT mix old with new batteries nor mix batteries of different compositions/brands.

-Do not heat or use deformed batteries.

-Do not dispose of batteries in normal waste.

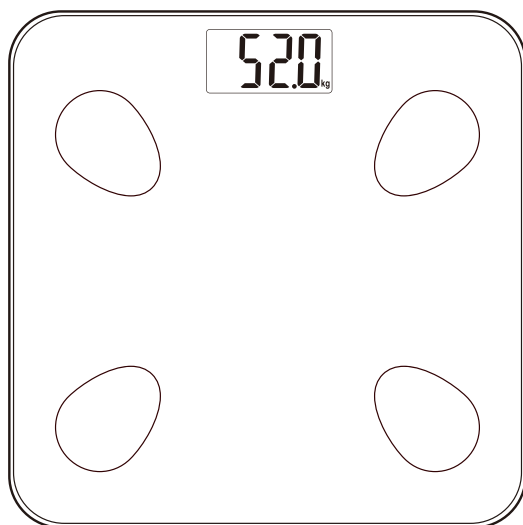-Check with your local authority for battery recycling advice.

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