
Water Flosser FP1

INSTRUCTION MANUAL

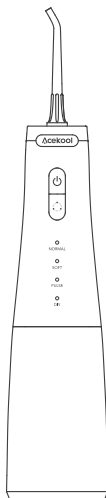EN

DE

FR

ES

IT

Thank you for purchasing Acekool Water Flosser. Before operating this unit, please read these instructions completely and keep the manual ready for further reference.

Contents

EN English 01-14

DE Deutsche 15-29

FR Française 30-44

ES Español 45-59

IT Italiana 60-74

EN Contents

■ Important Safeguards	02
■ FCC Statement	03
■ Safety Measures	04
■ Package Contents	06
■ Accessories	07
■ Specifications	07
■ Steps	08
■ Cleaning and Maintenance	12
■ Environmental Protection	13
■ Troubleshooting	14

■ ■ Important Safeguards

When using electrical products, especially when children are present, basic safety precautions should always be followed, including the following:

READ ALL INSTRUCTIONS BEFORE USING.

DANGERS

To reduce the risk of electrocution:

- Do not use while bathing.
- Do not place or store product where it can fall or be pulled into a tub, shower stall or sink.
- Do not place in or drop into water or other liquid.
- Do not reach for a product that has fallen into water. Unplug immediately.
- Do not handle charger with wet hands when charging.

WARNING

To reduce the risk of burns, electrocution, fire, or injury to persons.

- A product should never be left unattended when plugged in.
Exception: Contact lens disinfectors, toilet seat assemblies, and similar products.
- Close supervision is necessary when this product is used by, on, or near children or invalids.
- Use this product only for its intended use as described in this manual. Do not use attachments not recommended by the manufacturer.
- Never operate this product if it has a damaged cord or plug, if it is not working properly, if it has been dropped or damaged, or dropped into water. Return the product to a service center for examination and repair.
- Keep the cord away from heated surfaces.
- Never block the air openings of the product or place it on a soft surface, such as a bed or couch, where the air openings may be blocked. Keep the air openings free of lint, hair, and the like.
- Never use while sleeping or drowsy.
- Never drop or insert any object into any opening or hose.
- Unplug this product before filling. Fill (reservoir) with water only unless otherwise specified by manufacturer. Do not overfill (or specify filling instructions)

- The instructions for appliances having a part of class III construction supplied from a detachable power supply unit shall state that the appliance is only to be used with the power supply unit provided with the appliance.
- This product needs to be charged with the charger specified by the manufacturer, do not use the charger not specified or other charger.
- This appliance contains batteries that are non-replaceable.
- This appliance is not intended for use by persons (including children) with reduced physical, sensory or mental capabilities, or lack of experience and knowledge, unless they have been given supervision or instruction concerning use of the appliance by a person responsible for their safety.
- Children should be supervised to ensure that they do not play with the appliance.
- Use this product only for its intended use as described in this manual. Do not use attachments not recommended by the manufacturer.
- Do not use outdoors or operate where aerosol (spray) products are being used or where oxygen is being administered.

■ ■ ■ FCC Statement

This device complies with Part 15 of the FCC Rules. Operation is subject to the following two conditions: (1) this device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.

Warning: Changes or modifications to this unit not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment.

NOTE: This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to Part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio / TV technician for help.

Safety Measures

Basic safety precautions should be observed during the use of electrical products, especially when children are present.

Please read all safety guidelines before using this device.

WARNING

To reduce the risk of electric shock and other hazards, please observe the safety guidelines below:


- Check if the charger or power cord are damaged before initial and during continued use of this device. If either are damaged, do not attempt to use the device. Instead, contact after-sales services for a replacement.
- Do not throw the device into a fire or other heated environments.
- Do not immerse the device into water, as this may cause damage.
- Do not attempt to plug or unplug the charger while your hands are wet. When plugging or unplugging the device, do so using the charger, not the power cord.
- Do not modify, damage, or destroy the charging equipment. Do not forcefully stretch, bend, or twist the power cord. Do not place heavy objects on the power cord or stick the power cord between other objects.
- Do not drop or insert foreign objects into any of the device openings.
- Do not place the device, charger, or power cord on heated surfaces.
- Do not place the device in a hot or humid location, nor any locations exposed to direct sunlight.
- Before cleaning and maintenance, ensure that the device is turned off and not charging.

 WARNING


To reduce the risk of burns, electric shock, fires, and other injury hazards, please observe the safety guidelines below.

- Do not use this device in any manner other than those outlined in the manual.
- Do not use this device while a nozzle is not installed.
- Persons (including children) without the full mental capacity or relevant experience and knowledge must clean and maintain this device under the supervision of a guardian.
- Place this device outside of the reach of children and infants. Do not allow this device to be used as a toy.
- Users suffering from periodontitis or those who have undergone oral surgery within the past 2 months should consult a doctor before using this device.
- Only the nozzles and accessories included with, sold with, or recommended by this device should be used.
- Do not disassemble or swap out parts from this device.
- To prevent burns, do not fill the tank or attempt to clean your mouth with water over 40°C.
- Do not use this device for any purpose other than oral hygiene.
- To prevent injury, do not press the nozzle forcefully into your gums or teeth.
- Do not use this device in outdoor locations where aerosols (sprays) or oxygen supplies are in use.
- Remove all cosmetic teeth accessories before using this device.
- Do not use iodine, salt water, or concentrated essential oils that are insoluble in water with this device. Doing so will reduce its performance and shorten its service life.
- Do not fill the device with pure mouthwash liquid. If you wish to use a mouthwash, dilute it in accordance with the proper proportions.

■ ■ Package Contents

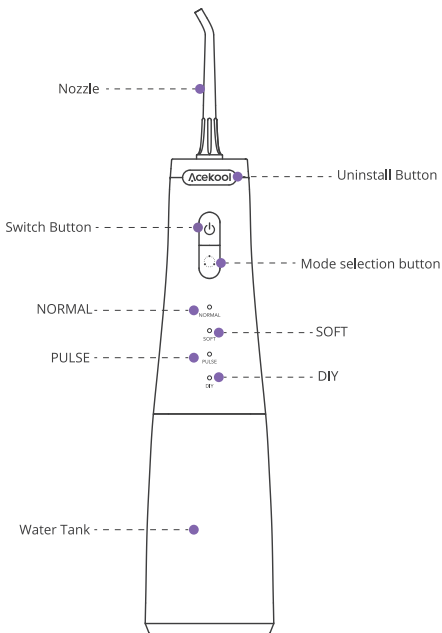Accessories

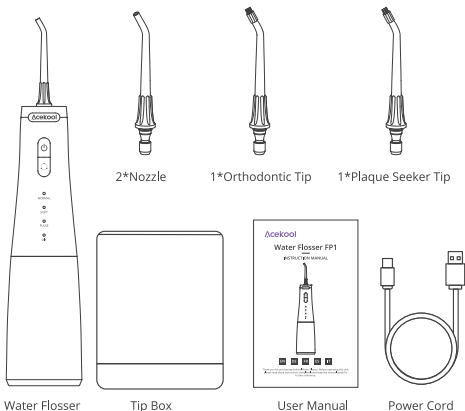Tip: Replace every 6 months Standard Nozzle

Specifications

Rated Voltage	3.7V ===
Power Rating	8W
Waterproof Rating	IPX7
Charge Time	≤4H
Input	5V===1A
Adaptor Rating Value Input	100-240V~50/60Hz 0.3A
Adaptor Rating ValueOutput	5V===1A 5.0W

Steps


Before Starting: Very Important! Ensure the device is **fully** charged before using it for the first time by connecting the USB cable to the charging port and turning on the power. This device cannot be used while charging.

Step 1: Fill the Water Tank

Open the lid and then fill water into the reservoir. If the lid falls off, you can insert the gasket of the lid to the opening, then close the lid.

Note: Ensure that the device is turned off before filling the water tank.

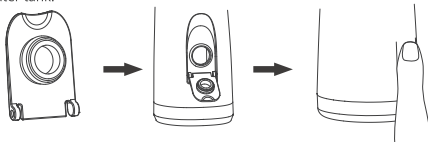Step 2: Install the Nozzle

Insert the nozzle into the hole on the top of the controller. You should hear a "click" sound on successful installation. Rotate the nozzle to adjust the direction of water outflow. If you need to replace the irrigator, press and hold the **uninstall** button, then remove the irrigator.

Installation

Replacement

Step 3: Select Mode



This device has 4 modes. Select the desired mode by tapping the mode button while the device is off. The functions of each mode are outlined below.

- After you press the on/off button, the device will enter the work mode. You can press the mode button to switch the mode in the following cycle: Normal level-Soft level-Pulse level-DIY. In any work mode, you can press and hold the on/off button for more than one second to turn off the device.

Enter the DIY mode method:

- Long-press the DIY mode button, and the DIY light flashes
- Water pressure circulates from the lowest to the highest, and then short-press the function conversion button at the appropriate water pressure
- Water pressure can be fixed.

- During charging, the function indicator displays the charging status, the battery displays the current power, and the charging indicators flash in sequence. When the device is fully charged, the battery power is full, and all the charging indicators are on. You can press and hold the mode button to check the current battery power, which is indicated by the number of indicators that are on. If the charging indicator flashes during use, you need to charge the device. When the device is nearly out of power, the charging indicator will flash for three times, then the device will be off automatically.
- The device will work for two minutes. Memory mode: When the device is on, it will memorize the last mode you use for more than one minute before you turn off the device.
- Hold the device upright during use to make sure that it is in the optimal working condition. If the device operates with any mode for more than 120 seconds, it will be off automatically.

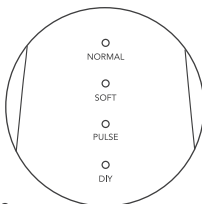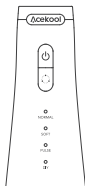Step 4: Before Use

Before using the device, lean slightly over the sink, extend the nozzle into your mouth and aim it at the teeth, gums, or between the teeth. Rotate the nozzle to adjust the direction of the water outflow so that it is perpendicular to the part of your mouth you wish to clean. Ensure that your mouth is open enough so that any water flows right out.

Note:



When using the oral irrigator for the first time, start from the lowest gear. Water pressure can be gradually increased as the oral cavity adapts. You may experience itching or minor gum bleeding during use. This is normal and should typically stop after 1-2 weeks.

Step 5: Use

Press and hold the switch button for 0.6 second to start the device. To pause the device during use, tap the switch button. Press and hold the button again to resume operation. This device is equipped with a 120-second shutdown timer. Automatic shutdown will only trigger if the water tank is empty. If you wish to keep using the device, tap the switch button to turn it off, then restart the device after refilling the water tank. Turn off the device after you are finished using it.


Techniques

For best results, start with the molars (posterior teeth) and start cleaning towards the front teeth. While cleaning, move the nozzle 90° along the gum line. Stay between each tooth for a short period. Close your lips slightly to avoid the water from spraying out, but leave enough space so that water can flow freely from your mouth to the sink.

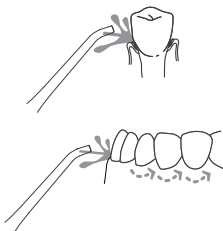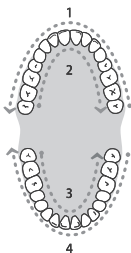Step 6: After Use

After use, open the water tank and pour out any remaining liquid. Turn on the device and drain any liquid left until the irrigator no longer sprays out any liquid, then turn off the device.

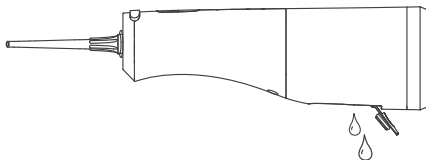Step 7: Charging

The indicator light is white, and will remain on while the device is in use. If the device shuts down after the indicator flashes for five times, the device is in low-voltage state. In the low-voltage state, after you press the on/off button again, the device will not power on. The indicator flashes for four times, indicating that you need to charge the device in time. When the device is charging, the indicator flashes. While the device is fully charged, the indicator light will remain on.


Cleaning and Maintenance

Cleaning the Body:

- Clean this device only with clean water or a neutral detergent. Do not use corrosive or abrasive detergents. Do not use water above 50°C to clean any part of this device, as this may lead to damage or reduced performance.
- Do not soak the oral irrigator in water. If necessary, clean the surface of the device by wiping with a soft cloth dipped in mild, non-corrosive detergent.

Cleaning the Tank:

- You can hold the device with your left hand, seize the buckle of the cover with your right hand, rotate the buckle left and right, then gently pull it down.
- Rinse the inside of the tank with clean room temperature water. Do not clean the tank with any type of detergent.
- After cleaning, locate the bottom cover to the buckle, then gently insert it to the device.

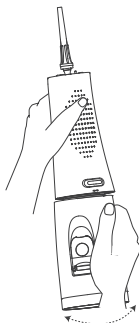Removing Hard Water:

- Depending on the mineral content of the water used, hard water deposits may build up on the device. Failure to remove hard water buildups will reduce product performance.
- Point the nozzle at the sink, turn on the device, then let it run until all water has been used. Fill the water tank with clean warm water, then repeat the above steps to clean the tank. This cleaning process should be repeated every 1-3 months to maintain product performance.


Environmental Protection

Battery Maintenance

- Ensure that the battery is fully charged before initial use.
- Ensure that the battery is fully charged before storing the device for an extended period.
- In order to prolong battery life, charge the device immediately when the battery is low.

Battery Recycling

This device uses recyclable rechargeable batteries. At the end of the device's lifecycle, recycle or dispose of the batteries separately in accordance with your local waste management regulations before disposing of the device.

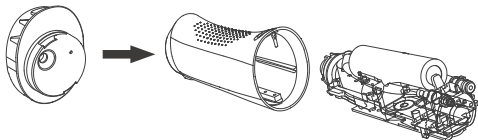Recycling Procedure: While the device is turned off, remove the nozzle and water tank, then use a screwdriver or a similar tool to pry open the base to remove the recyclable batteries.

Note: Do not short circuit the positive and negative terminals of the removed battery. The terminals can be insulated with tape.

■ ■ Troubleshooting

Problem	Cause	Solution
There is a leak between the nozzle and the body	The seal of the nozzle holder is damaged, or the nozzle was not installed in the proper position	Contact seller for replacement, or reinstall the nozzle in the proper position
There is a leak in the water inlet	The seal is damaged, or there are impurities on the inlet	Contact seller for replacement
There is a leak in the water tank	The seal is damaged	Replace the seal, or contact seller for replacement
There is not enough water pressure	Not enough charge	Try again after charging the device
The device cannot start	Not enough charge	Try again after charging the device, 12 or contact seller for replacement
The device is not releasing any water	The device is being used for the first time, or has not been used in an extended period	Add a suitable amount of water to the nozzle installation hole and try again after reinstalling the nozzle.

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Distributor: Acekool Innovation Inc.