

尺寸：170mm\*100mm

骑马钉

封面128g铜版纸，内页80g书纸


ANJOU | [www.anjou.com](http://www.anjou.com)  
FIND YOUR BEAUTY

ANJOU

# Beauty Machine

AJ-PCA050

User Manual

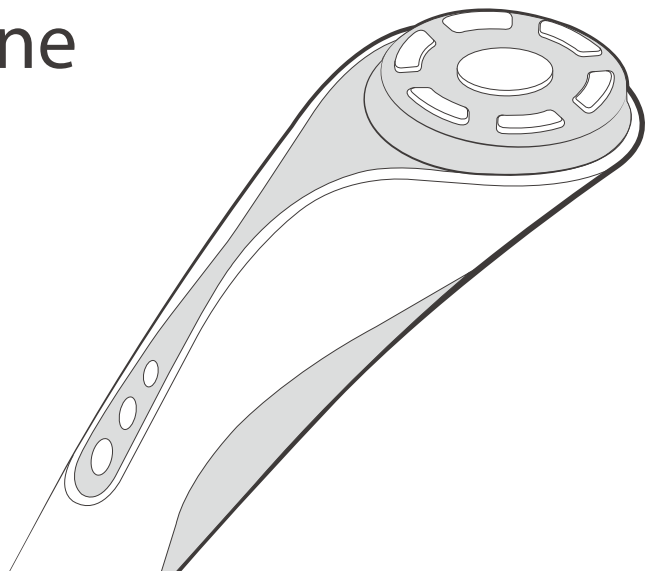

# CONTENTS

---

01-02 SAFETY INSTRUCTIONS

03-04 KNOWING YOUR BEAUTY MACHINE

- 03 Beauty Machine
- 04 Accessories
- 04 Specifications

05-13 USING YOUR BEAUTY MACHINE

- 05 Skin Allergy Test
- 06 Treatment Modes & Frequency
- 07 Intensity Levels
- 07 Light Treatment Effects
- 08 Treatment Skills
- 10 Daily Use


14 CLEANING AND MAINTENANCE

- 14 Cleaning
- 14 Store & Disposal

15 TROUBLESHOOTING

16 WARRANTY

# SAFETY INSTRUCTIONS

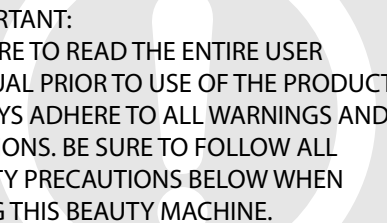

**IMPORTANT:**  
BE SURE TO READ THE ENTIRE USER  
MANUAL PRIOR TO USE OF THE PRODUCT.  
ALWAYS ADHERE TO ALL WARNINGS AND  
CAUTIONS. BE SURE TO FOLLOW ALL  
SAFETY PRECAUTIONS BELOW WHEN  
USING THIS BEAUTY MACHINE.

## Related to Health

1. This instrument is ONLY used for healthy skin. DO NOT use on plastic surgical spots, wounds or ulcers.
2. If your skin is sensitive, please perform a skin allergy test before using the instrument.
3. DO NOT use during periods of unstable hormone secretion such as menstrual period, pregnancy, and lactation.
4. Those who have a pacemaker, a history of epilepsy or heart disease, metal implants, or other serious diseases should not use the instrument.

## Electricity Safety

1. ONLY use the original adapter to power up the instrument.
2. ALWAYS remove metal accessories from your body before use if any.
3. Be sure that your hands are completely dry when operating the instrument.
4. Keep away from sink or bathtub when using. DO NOT submerge the instrument under water.
5. DO NOT touch the instrument should it fall into water. Unplug the instrument immediately with dry hands.
6. Turn the instrument OFF and then unplug the adapter from the socket after use or before cleaning and maintenance.

## Product Safety


1. DO NOT disassemble the instrument by yourself.
2. Keep away from places where children can reach the instrument.
3. Keep away from fire or flammable materials.
4. Avoid falling from height or any impact.

## Related to Use


1. NEVER let the instrument stay on skin too long during treatment to avoid overheat.
2. If there is any discomfort during use, stop using immediately and consult a professional medical practitioner.
3. If the instrument is damaged or does not work normally, please contact Anjou Customer Support.
4. The instrument is not a medical instrument, and CANNOT be managed and used for medical use.

# KNOWING YOUR BEAUTY MACHINE

## Beauty Machine

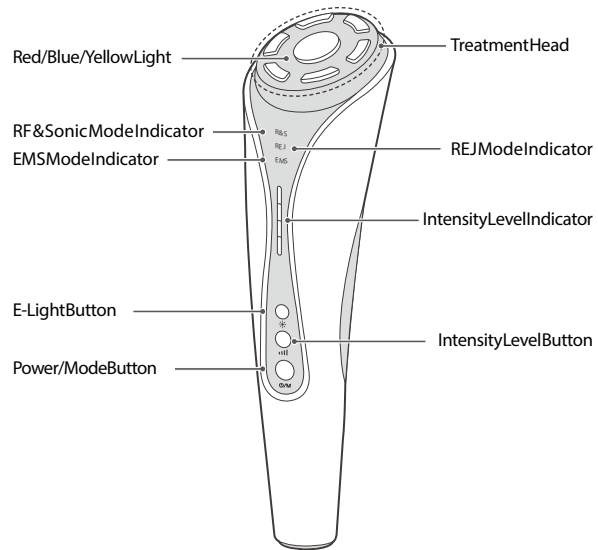

## Accessories

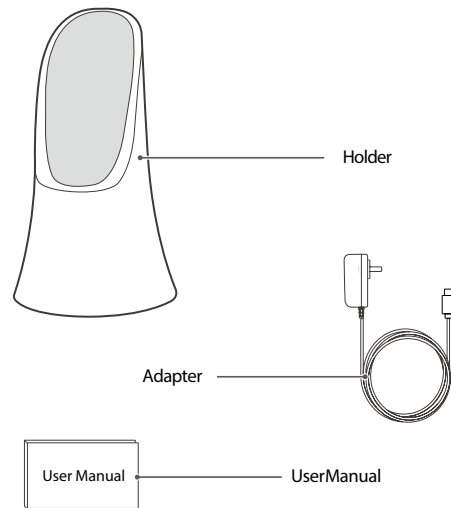


## Specifications

Model	AJ-PCA050
Input	100-240V
Output	9V = 1500mA

# USING YOUR BEAUTY MACHINE

## Skin Allergy Test

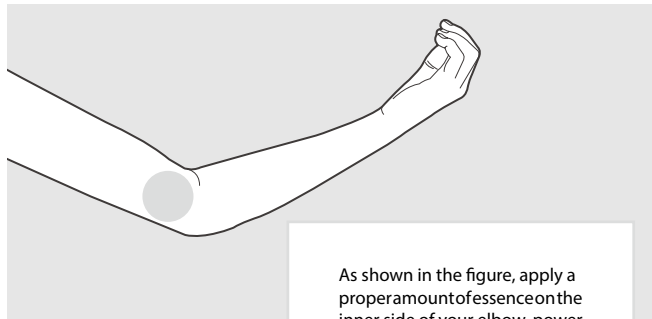

As shown in the figure, apply a proper amount of essence on the inner side of your elbow, power on the instrument, select a mode and a proper intensity level, massage the area for 2-3 mins, power off the instrument, and air dry the area. If no allergic reaction within 48 hours, you can use this instrument normally.

## Treatment Modes & Frequency

**EMS Mode (Tighten)**  
Hold the **⏻/M** button for 1.5s (Power on) to enter EMS mode.



Default Light: Red & blue



Mode Effects: Micro-current exercises muscles, eliminates edema of the skin, lifts the jawline, tightens the skin

**REJ Mode (Deep SPA)**  
After powering on, press the **⏻/M** button once to enter REJ mode.



Default Light: Red & yellow



Mode Effects: Innovative Multipolar RF, brand new SPA experience, deeply rejuvenates your skin, promotes collagen production, anti-aging, and wipes out age signs

**R&S Mode (RF & Sonic)**  
After powering on, press the **⏻/M** button twice to enter R&S mode.



Default Light: Red



Mode Effects: Introduces nutrition from deep inside to skin surface, enhances skin elasticity and stimulates collagen production

#### Treatment Frequency

Use the instrument according to the following frequency, better results can be obtained. You can reduce the frequency after two weeks.

Mode	Week 1-2	Week 3-8	After Week 8
EMS	Once a day 4-8 mins	Once every 2 days 4 mins	Once a week 4 mins
REJ	Once every 2 days 4-8 mins	Once every 2 days 4 mins	Once a week 4 mins
R&S	Once every 2 days 4-8 mins	Once every 2 days 4 mins	Once a week 4 mins

## Intensity Levels

There are 4 levels for each mode. During your first treatment session, we recommend that you start with level 1 and then adjust the level according to your skin's tolerance for subsequent treatments.

## Light Treatment Effects

#### Red Light Treatment Effects:

The wavelength of red light is  $622\pm 5\text{nm}$ . It can go deep into the skin, activate fibroblast cells, accelerate the replacement of cells, and stimulate collagen regeneration. Long-term using of red light can smooth wrinkles and fine lines, effectively reduce eyebags, and make skin smooth and elastic.

#### Blue Light Treatment Effects:

The wavelength of blue light is  $467\pm 5\text{nm}$ . It can eliminate the microorganism that causes acne, and relieve inflammation and skin redness. Long-term using of blue light can improve the water-oil balance of your skin, reduce acne and other skin issues, and keep skin in a healthy status.

#### Yellow Light Treatment Effects:

The wavelength of yellow light is  $580\pm 3\text{nm}$ . It is mild so that can be used for sensitive skin. It can enhance skin immunity by moderately stimulating the lymphatic system and the nervous system. Long-term using of yellow light can improve dull skin, fade spots, brighten and whiten your skin.


- We suggest you to treat with R&S mode first to introduce nutrition before treatment with other modes.
- Each mode is well-designed with its specific default LED color for the best treatment effect, however, you can change or turn off the LED color light as per your preference.

## Treatment Skills


#### Tips:

- Massage the area marked with numbers in order and stay 3 seconds for each one.
- NEVER let the instrument stay on skin too long during treatment to avoid overheat.


A. Mouth corners drooping, nasolabial folds obvious (use REJ - deep SPA mode)

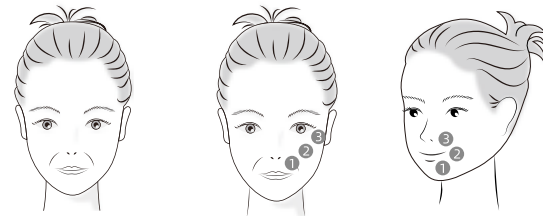


B. Tighten your face (use REJ. - deep SPA mode)

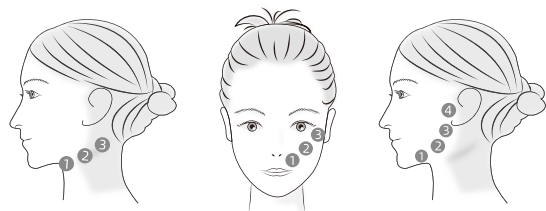

C. Periorbital and forehead wrinkles (use R&S – RF ultrasonic treatment mode)

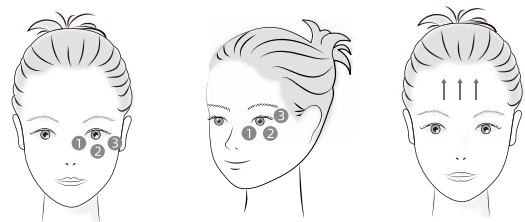


D. Small facial muscle exercises (use EMS - tighten mode)

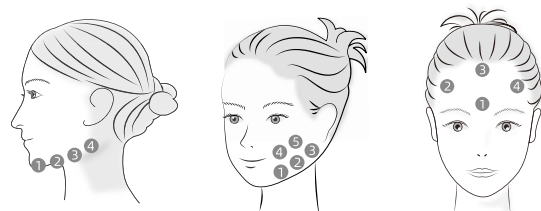


Daily Use (Note: Essence is not included in the package.)

1




Clean the skin area you want to treat, and evenly apply a proper amount of essence on the skin.

2

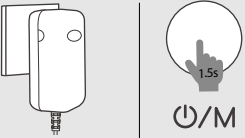

Plug the adapter into a socket and hold the **⏻/M** button for 1.5s to power on.

3

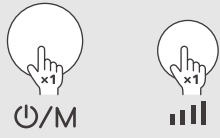


Press the **⏻/M** button to select the desired mode. Press the **||||** button to select the desired intensity level.

4



Keep the treatment head cling to skin to ensure proper heating of the treated area. Slowly and constantly move the instrument over the skin.

5

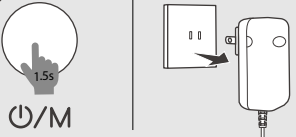

Hold the **⏻/M** button for 1.5s to power off and unplug the adapter.


6




Wipe down the treatment head with a damp cloth or paper towel.




Regular Treatments (face)

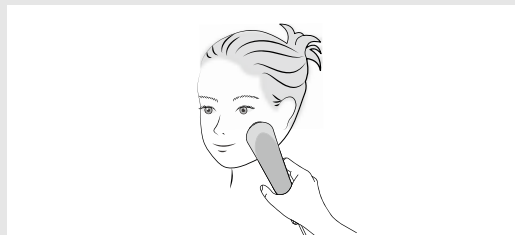

From undernose to eyebrow corner

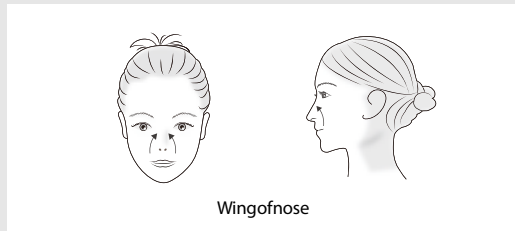


Wing of nose

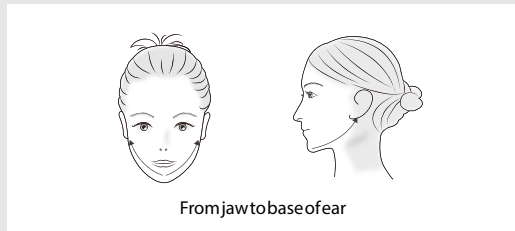


From jaw to base of ear


Regular Treatments (face)

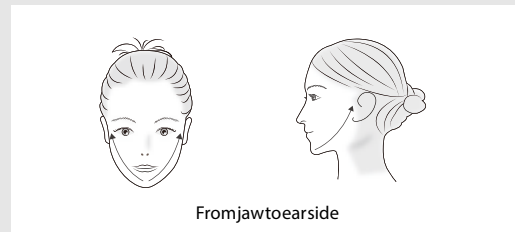


From jaw to ear side

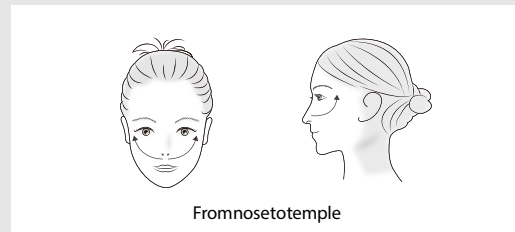

From nose to temple

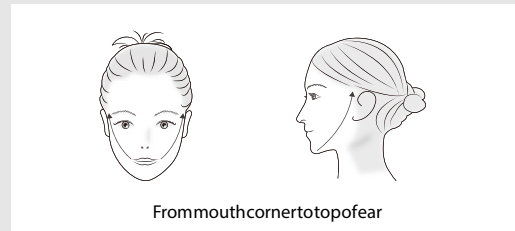


From mouth corner to top of ear

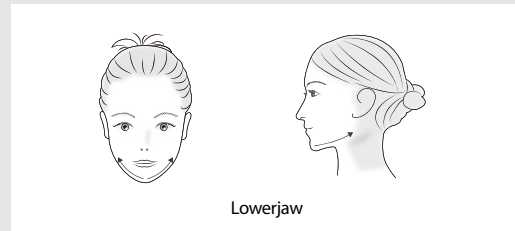


Lower jaw

### Regular Treatments (neck and others)

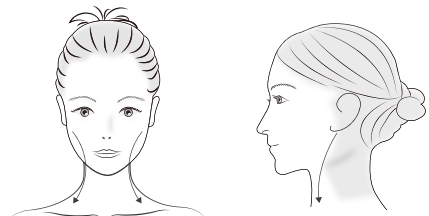

From base of ear to collarbone

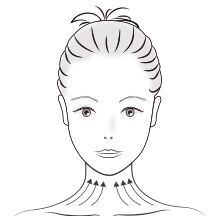

From top of collarbone to face line (as shown in figure).  
Time: 2 mins

You can treat arms, waist & lower abdomen, foot & calf, hip and hands as well.



#### Note:

- DO NOT rinse the treatment head or the instrument with water.
- The instrument is powered up by AC only and non-rechargeable.

## CLEANING AND MAINTENANCE

### Cleaning

1. After treatment, power off and unplug the instrument.
2. Wipe down residue on the treatment head with a damp cloth or paper towel.
3. Place the instrument back in the holder.



DO NOT rinse the instrument under tap water.



DO NOT clean with any chemical detergents.



You CANNOT disassemble the instrument by yourself.

### Store & Disposal

Store the instrument in a dry and shady place.  
DO NOT dispose this product as general trash. Follow your local government rules for proper disposal of this product.

## TROUBLE-SHOOTING

Q: Will I get acne if I use essence with the instrument? Do I need to clean after use?

A: It will not cause acne if the essence used is suitable for your skin. To prevent from excess nutrition and your skin doesn't absorb, you'd better clean the skin with water after use, and then moisturize.

Note: Essence is not included in the package.

Q: Is the using time of each mode extendable?

A: To avoid burdening the skin, each mode of the instrument is set with an intelligent timer. It is not recommended to extend the using time.

Q: Can I use the instrument during menstrual period?

A: Since the secretion is unstable during the menstrual period, it is not recommended to use during the menstrual period. Use the instrument after the menstrual period.

Q: Can I proceed to the next treatment after one treatment is over?

A: Yes. There is no need for any recovery period or rest time if you have no discomfort nor pain during use.


Q: My skin is sensitive. Will it be allergic or painful?

A: If so, you can do the allergy test before use. You might feel a slight pain when using EMS mode, which is normal. It is recommended to use the lowest intensity for sensitive skin. Slight redness or heating after use is normal. If adverse reactions such as allergies occur, please stop using it immediately.

## WARRANTY

Thank you for choosing Anjou Beauty Machine. This product is covered with Anjou product and labor warranty for 12 months from the date of its original purchase.

We can only provide after-sales service for products that are sold by Anjou or Anjou authorized retailers and distributors. If you have purchased your unit from a different place, please contact your seller for return and warranty issues.

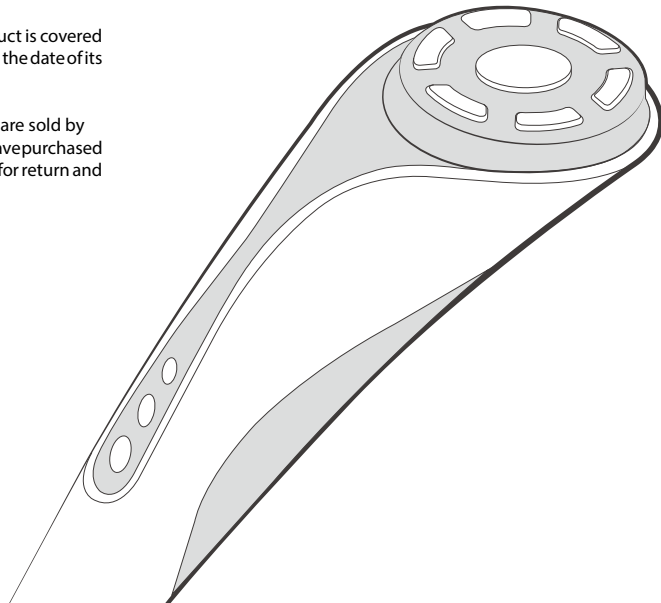



## WEEE Compliance

Correct Disposal of This Product (Waste Electrical & Electronic Equipment)

(Applicable in countries with separate collection systems)

This marking on the product, accessories or literature indicates that the product and its electronic accessories should not be disposed of with other household waste at the end of their working life. To prevent possible harm to the environment or human health from uncontrolled waste disposal, please separate these items from other types of waste and recycle them responsibly to promote the sustainable reuse of material resources. Household users should contact either the retailer where they purchased this product, or their local government office, for details of where and how they can take these items for environmentally safe recycling. Business users should contact their supplier and check the terms and conditions of the purchase contract. This product and its electronic accessories should not be mixed with other commercial wastes for disposal.

## FCC Compliance

This device complies with Part 15 of the FCC Rules. Operation is subject to the following two conditions:

(1) This device may not cause harmful interference, and (2) This device must accept any interference received, including interference that may cause undesired operation.

 [support@anjou.com](mailto:support@anjou.com)

 [www.anjou.com](http://www.anjou.com)



Manufacturer: Shenzhen Danya Technology Co.,Ltd.  
Address: Room 705, Building E, Galaxy World Phase II, Minle  
Community, MinZhiStreet, LonghuaDistrict, Shenzhen, China

MADEINCHINA

