

breo®

breo®

iSee 4

EYE MASSAGER



USER MANUAL

Thank you for purchasing BREO iSee4.
Please read this manual carefully before using the product,
and keep it in a safe place for future reference.


Contents

Inside the Package	2
Function Deployment	3
Preparation	4
Operation	5
Precautions	6-8
Troubleshooting & Technical Parameters	9

Inside the Package

The package includes: 1 Main Unit, 1 Power Adaptor, 1 USB Cable, 1 Travel Case, 1 User Manual

Massager

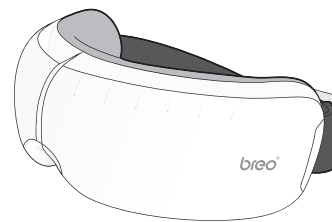

Main Unit


Basic components



1 USB Cable 1 Power Adaptor User Manual Travel Case

Function Deployment

1. Power Indicator
2. Main Button
3. LCD Screen
4. Adjustable Elastic Band
5. DC Power Port

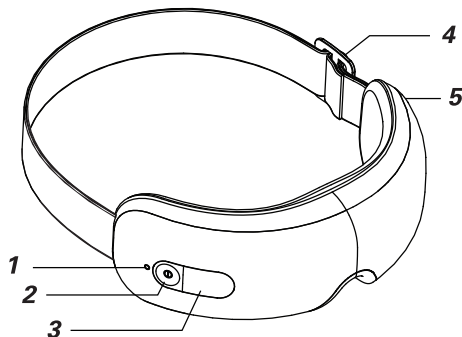

3


Preparation

Prior to use, check the condition of the machine and accessories:

- a. No damage on power charger
- b. Battery power is sufficient


Charging

1. Open the cover of the DC power port on the main unit.
2. Insert the plug of the power charger to the DC port.
3. During charging, the battery icon on the LCD screen flashes and shows the increasing battery power. The icon stops flashing when fully charged.
4. The massager cannot be used while charging. If the unit is turned on while being charged, it will prompt the self-protection program and the massager will automatically turn off.
5. The recommended charging duration is 2-3 hours.

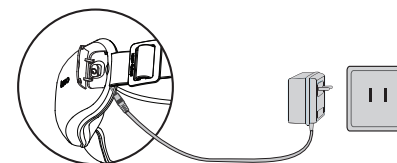

Wear

1. Put on the main unit where it will cover the eyes.
2. Adjust the elastic band to achieve the most comfortable fit.

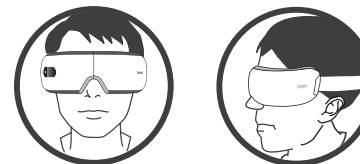

4


Operation

Power On

Press and hold the ① Button.
Preset default is SLEEP mode.

Mode Selection

Quick Press the ① Button to select mode.

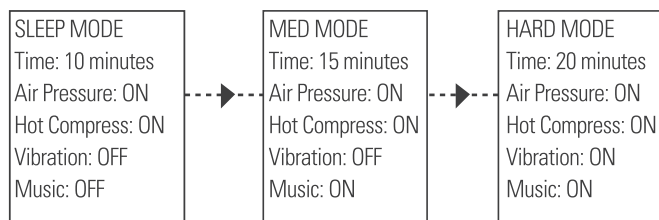

Power Off

Press and hold the ① Button again

Precautions

Please read this section carefully for correct and safe use of this product.

1. Please read the manual carefully and operate accordingly.
2. Children under 7 years old cannot use the massager.
3. This massager should not be used by:
 - a. People who suffer from brain disease or have undergone brain surgery
 - b. People who are suffering from glaucoma, retinoschisis or other serious eye diseases
 - c. People who have physical or mental disorder
4. Take off contact lenses/eye glasses before use.
5. Put on/wear the unit first BEFORE turning the power on.
6. Do not use while walking or driving.
7. Keep your eyes closed for a few seconds after use, then open them slowly.
8. Keep the massager away from humidity, heat, and/or corrosive gas.
9. If the massager will not be used for a long time, remove the batteries to avoid battery leakage.
10. This massager features heat compression. Individuals who are incapable of or do not ordinarily feel heat/cold sensation should use this device with caution.
11. To clean, wipe with aseptic cotton and alcohol. Do not wash or immerse in water.
12. This product is not a substitute for healthcare professionals' advice and treatment. In case of doubt or if necessary, please consult your healthcare professional.

Precautions

13. This device complies with Part 15 of the FCC Rules. Operation is subject to the following two conditions:
 - (1) this device may not cause harmful interference, and
 - (2) this device must accept any interference received, including interference that may cause undesired operation.
14. Please note that changes or modifications not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment.
15. This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to Part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:
 - Reorient or relocate the receiving antenna.
 - Increase the separation between the equipment and receiver.
 - Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
 - Consult the dealer or an experienced radio/TV technician for help.

Precautions

16. **IMPORTANT:** Any individual who may: be pregnant or recently gave birth, has a heart condition or a pacemaker, suffers from diabetes or other sensory impairment, phlebitis and/or thrombosis, joint dysfunction, inflammation or injuries, weak bones, osteoporosis, an abnormal or curved spine, is at an increased risk of developing blood clots, has pins/screws/artificial joints or other medical devices implanted in his/her body, is being treated for a medical condition or physical injury or is under a doctor's care should consult with a medical professional before using a massaging device designed for home use.


Troubleshooting & Technical Parameters

When operational problems occur, please refer to the chart below.

Do not disassemble the machine or try to repair it by yourself.

Problem	Possible Reason	Possible Solution
Power-ON failure	Lithium battery power is low	Recharge it at once
Weak air pressure or vibration	Lithium battery power is low	Recharge it at once

Specifications

Power Supply	DC: 5V  1A
Power	5W
Ideal Environment	Temperature: 5°C-40°C Relative Humidity: ≤80% Atmospheric Pressure: 860hpa ~ 1060hpa
Dimension	Main Unit L 8in x W 3in x H 4in
N.W.	11.2 OZ

breo®

Made in China

Manufactured by: Body Source Technology (ShenZhen) Co.,LTD
Company Address: 8thFloor, Building H, Gang Zhi Long Technology Park,
He Ping East Road, Long Hua Street, BaoAn, ShenZhen, 518000, CHINA