




Visual Water Flosser

- User Manual -

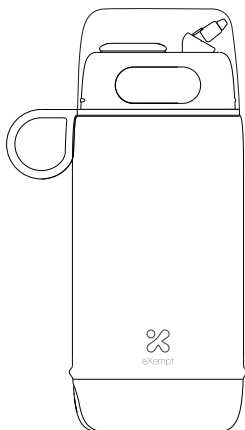

Model: C3 Pro

Important Safeguards



Read All Instructions before Using

When using electrical products, especially when children are present, basic safety precautions should always be followed, including the following:

To reduce the risk of electrocution:

- Do not handle charger with wet hands.
- Do not place charger in or drop charger into water or other liquid.
- Check the charger cord for damage before the first use and during the life of the product.

To reduce the risk of burns, electrocution, fire, or injury to persons:

- Do not plug this device into a voltage system that is different from the voltage system specified on the device or charger.
- Do not charge this unit with any charging devices other than the one provided.
- Do not use this product if it has a damaged cord or plug, if it is not working properly, if it has been dropped or damaged, or dropped into water.
- Do not direct water under the tongue, into the ear, nose or other delicate areas. This product is capable of producing pressures that may cause serious damage in these areas.
- Do not direct water into the nose or ear. The potentially deadly amoeba, *Naegleria fowleri*, may be present in some tap water or unchlorinated well water and may be fatal if directed into these areas.
- Use this product only as indicated in these instructions or as recommended by your dental professional.


- Only use tip heads and accessories that are recommended by eXempt.
- Do not drop or insert any foreign object into any opening on the product.
- Keep charging unit and cord away from heated surfaces.
- Do not use iodine, saline (salt solution), or water insoluble concentrated essential oils in this product. Use of these can reduce product performance and will shorten the life of the product. Use of these products may void warranty.
- Remove any oral jewelry prior to use of this product
- Do not use if you have an open wound on your tongue or in your mouth.
- If your physician advised you to receive antibiotic premedication before dental procedures, you should consult your dentist before using this instrument or any other oral hygiene aid.
- This appliance is not intended for use by persons [including children] with reduced physical, sensory or mental capabilities, or lack of experience and knowledge, unless they have been given supervision or instruction concerning use of the appliance by a person responsible for their safety
- Children should be supervised to ensure that they do not play with the appliance
- If the supply cord is damaged, it must be replaced by a special cord or assembly available from the manufacturer or its service agent.

Product Introduction

Congratulations on your new eXempt visual water flosser! With 500W HD endoscope technology and 3 mode settings, the flosser can let you personalize your water pressure intensity with soft, strong and pulse options by seeing your teeth conditions on your phone. Plus it's collapsible design and tip heads storage which bring you unprecedented convenience.

Packing List


- eXempt Visual Water Flosser C3 Pro x 1
- Classic Tip Head x 1
- Pocket Tip Head x 1
- Orthodontic Tip Head x 1
- USB Cable x 1
- Quick Guide x 1

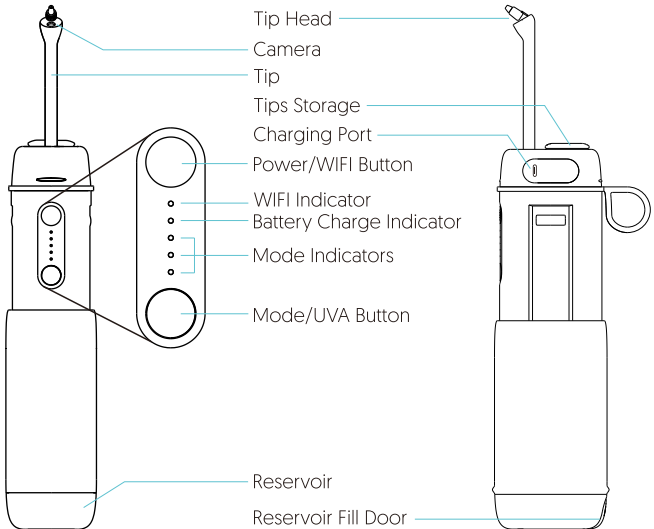


Features:

- 500W HD endoscope- provide clear and stable images of your teeth
- Collapsible design and tip heads storage- convenient for travel
- Replaceable three tip heads- classic tip head, pocket tip head and orthodontic tip head
- IPX7 waterproof rating- for better clean your flosser.
- USB-C charging- provide more power and quick charging


Using your eXempt Visual Water Flosser:

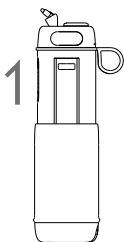


1. Pull out reservoir from base to the usage line.

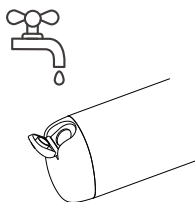



2. Lift open reservoir door and fill with lukewarm water.

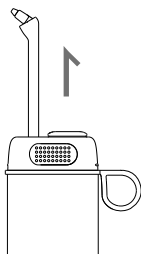

3. Pull out the tip head..

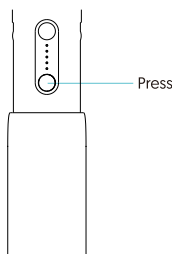

4. Press the mode button to choose a mode.

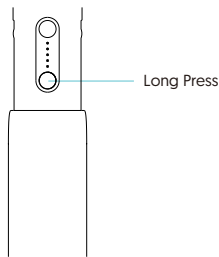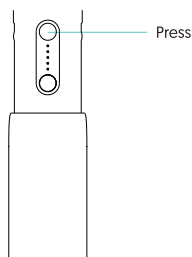

5. Press the power button to start flossing.

6. Long press the UVA button to turn on/off the UVA.

Mode Guide



- ① **Soft Mode:** standard for superior cleaning, suitable for sensitive teeth and beginners
- ② **Strong Mode:** ideal mode for removing plaque and debris deep between teeth and below the gumline.
- ③ **Pulse Mode:** gently massage your gums.

eXemptcares APP- Getting Started



1. Download the eXemptcares app to your phone.




eXemptCares





2. Open the app, tap “Add device” or “+”, Choose the visual water flosser to add the device.

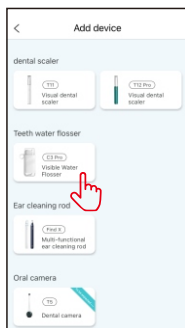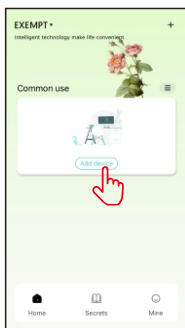

3. Click “Bind now”, then long press the power button to turn on WIFI, Next.

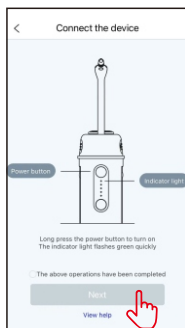

4. Click “start connecting” to enter into WIFI setting page of your phone, Pair your flosser by connecting WIFI “eXempt-C3 Pro-xxx”, then return to app after successfully connecting..

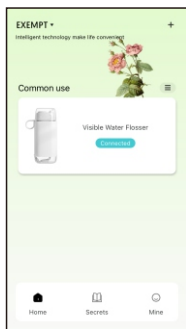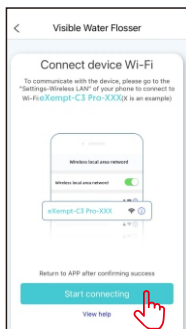

APP Icon Description



- **Undirect:** Flip the screen up and down
- **Static Mode:** Standard screen
- **Split Screen:** Double screen to display your teeth by the flosser otoscope
- ▮|▮ Flip the screen left and right
- Start flossing
- Shoot a photo or video
- ▮| Change your screen to noir display
- 🖼 Check your photos or videos
- Select photo and video shooting
- Select Modes

Warranty and Support

eXempt warrants its products for 2 years after the date of purchase. Defects due to faulty materials and workmanship will be replaced at Exempt's expense provided that convincing proof of purchase in the qualifying period is provided. Use of unauthorized replacement parts will void this warranty.

Warranty Exclusions

What is not covered under warranty

- Damage caused by use of unauthorized replacement parts or unauthorized flosser tips.
- Damage caused by misuse, abuse, neglect, alterations or unauthorized repair.
- Normal wear and tear, including chips, scratches, abrasions, discoloration or fading.

Contact Information

www.exemptcares.com

info@exemptcares.com

Search eXempt Cares

