

## Blood Pressure Monitor

Arm Type



**Healthcare-Manager.com**

Manufactured for Easy Healthcare Corporation

360 Shore Dr. Unit B

Burr Ridge, Illinois, 60527, UNITED STATES

Questions or comments? Please call toll-free:

1-855-822-6999 M-F 9 a.m.-5 p.m. CST

E-mail: [service@healthcare-manager.com](mailto:service@healthcare-manager.com)

**EBP-095**

**User Manual**

Table of Contents

INTRODUCTION..... 2

- General Description
- Safety Information
- LCD Display Signal
- Monitor Components

BEFORE YOU START..... 6

- The Choice of Power Supply
- Installing and Replacing the Batteries
- Setting Date, Time and Measurement Unit

MEASUREMENT..... 9

- Tie the Cuff
- Start the Measurement

DATA MANAGEMENT..... 11

- Recall the Records
- Delete the Records

INFORMATION FOR USER..... 14

- Tips for measurement
- Maintenances

ABOUT BLOOD PRESSURE..... 16

- What are systolic pressure and diastolic pressure?
- What is the standard blood pressure classification?
- Irregular Heartbeat Detector
- Why does my blood pressure fluctuate throughout the day?
- Why do I get a different blood pressure at home compared to the hospital?
- Is the result the same if measuring on the right arm?

TROUBLESHOOTING.....18

SPECIFICATIONS.....19

AUTHORIZED COMPONENTS .....20

CONTACT INFORMATION.....20

COMPLIED STANDARDS LIST.....21

FCC STATEMENT.....23

EMC GUIDANCE.....23

## ▼ General Description

Thank you for selecting **easy@home**® arm type blood pressure Monitor (EBP-095). The monitor features blood pressure measurement, pulse rate measurement and the result storage. The design provides you with two years of reliable service.

Readings taken by the EBP-095 are equivalent to those obtained by a trained observer using the cuff and stethoscope auscultation method.

This manual contains important safety and care information, and provides step by step instructions for using the product.

**Read the manual thoroughly before using the product.**

Features:

- 60mm\*80mm Digital LCD display
- Maximum 60 records
- Measuring during inflation technology








## ▼ Indications for Use

The Blood Pressure Monitor is digital monitors intended for use in measuring blood pressure and heartbeat rate with arm circumference ranging from 22 cm to 32 cm ( about 8¾"-12½" ) or 22cm to 42cm(about 8¾"-16½").

It is intended for adult indoor use only.

## ▼ Safety Information

The signs below might be in the user manual, labeling or other component. They are the requirement of standard and using.

	Symbol for "THE OPERATION GUIDE MUST BE READ"		Symbol for "TYPE BF APPLIED PARTS"
	Symbol for "DIRECT CURRENT"		Symbol for "ENVIRONMENT PROTECTION – Electrical waste products should not be disposed of with household waste. Please follow local guidelines."
	Caution: These notes must be observed to prevent any damage to the device		Symbol for "MANUFACTURE DATE"
	Symbol for "SERIAL NUMBER"		

## ⚠ CAUTION

This device is intended for adult use only.

This device is intended for non-invasive measuring and monitoring of arterial blood pressure. It is not intended for use on extremities other than the wrist or for functions other than obtaining a blood pressure measurement.

Do not confuse self-monitoring with self-diagnosis. This unit allows you to monitor your blood pressure. Do not begin or end medical treatment without asking a physician for treatment advice.

If you are taking medication, consult your physician to determine the most appropriate time to measure your blood pressure. Never change a prescribed medication without consulting your physician.

If the cuff pressure reaches 40 kPa (300 mmHg), the unit will automatically deflate. Should the cuff not deflate when pressures reaches 40 kPa (300 mmHg), detach the cuff from the arm and press the START/STOP button to stop inflation.

To avoid measurement errors, carefully read this manual before using the product.

The equipment is not AP/APG equipment and not suitable for use in the presence of a flammable anesthetic mixture with air with oxygen or nitrous oxide.

The operator shall not touch output of batteries/AC adapter and the patient simultaneously.

Do not wind air tube in the neck.

Please use ACCESSORIES and detachable partes specified/ authorised by MANUFACTURE. Otherwise, it may cause damage to the unit or danger to the user/patients.

Please note that Luer lock connectors are not used on the product and please DO NOT change any provided connectors.

Manufacturer will make available on request circuit diagrams, component parts list etc.

WARNING: No modifications of this equipment is allowed.

This unit is not suitable for continuous monitoring during medical emergencies or operations.

Otherwise, the patient's arm and fingers will become anaesthetic, swollen, and even purple due to a lack of blood.

Please use the device under the environment which was provided in the user manual. Otherwise, the performance and lifetime of the device will be impacted and reduced.

During use, the patient will be in contact with the cuff. The materials of the cuff have been tested and found to comply with requirements of ISO 10993-5:2009 and ISO 10993-10:2010. It will not cause any potential allergic reaction or contact injury.

The device doesn't need to be calibrated within the two years of reliable service.

Please dispose of ACCESSORIES, detachable parts, and the ME EQUIPMENT according to the local guidelines.

**When the device was used to measure patients who have common arrhythmias such as atrial or ventricular premature beats or atrial fibrillation, the best result may occur with deviation.**


This device may provide contradictory results for any female subject who may be suspected of, or is pregnant. Besides providing inaccurate readings, the effects of this device on the fetus are unknown.

The device has been evaluated clinically using manual cuff/stethoscope auscultation as the reference. Blood pressure measurements determined with this device are equivalent to those obtained by a trained observer using the cuff/stethoscope auscultatory method, within the limits prescribed by the American National Standard, Manual, electronic, or automated sphygmomanometers."

If you have any problems with this device, such as setting up, maintaining or using, please contact the SERVICE PERSONNEL of **easy@home**®. Don't open or repair the device by yourself.


Please report to **easy@home**® if any unexpected operation or events occur.

## ♥ LCD Display Signal

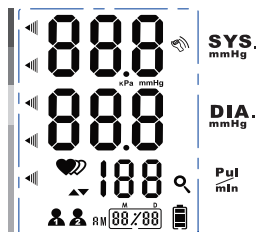

SYMBOL	DESCRIPTION	EXPLANATION
<b>SYS</b>	Systolic blood pressure	High pressure result
<b>DIA</b>	Diastolic blood pressure	Low pressure result
<b>Pul/min</b>	Pulse	Pulse/minute
	Deflating	CUFF air is exhausting of deflating
<b>RM 88/88</b>	Current Time	Time(year:month:day:hour:minute)
<b>M 88/88</b>	Memory	If "MEM" shows, the displayed measurement values is from the memory.
<b>kPa mmHg</b>	Measurement unit	Measurement Unit of the blood pressure (1mmHg=0.133kPa) (1kPa=7.5mmHg)
<b>Lo + </b>	Low battery	Batteries are low and need to be replaced
	Irregular heartbeat	Irregular heartbeat detection
	Grade	The grade of the blood pressure
	Heartbeat	Heartbeat detection during measurement
	Shocking reminder	Shocking will result in inaccurate
	User 1	Start measurement for user 1 and save the measuring result automatically
	User 2	Start measurement for user 2 and save the measuring result automatically
	Data Enquiry Mode	Recall the records

## ♥ Monitor Components

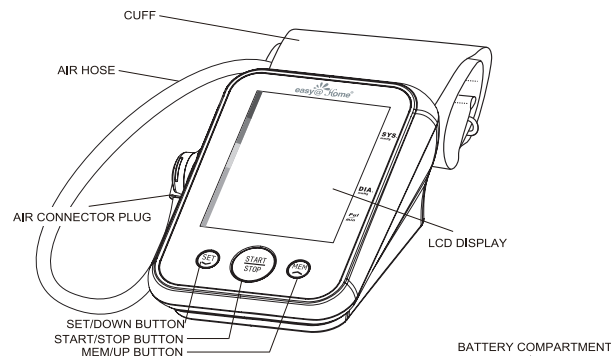

Component list of pressure measuring system

- 1 Cuff
- 2 Air pipe
- 3 PCBA
- 4 Pump
- 5 Valve

## ♥ List

1. Blood Pressure Monitor (EBP-095)



3. 4\*AA Batteries



- 2.Cuff (Type BF applied part)


(Please use **easy@home** authorized cuff 22cm~32cm (8"-12") or 22cm~42cm (8"-16"). The size of the actual cuff please refer to the label on the attached cuff.)



4. User manual

## ♥ The Choice of Power Supply

1. Battery powered mode:  
6VDC 4\*AA batteries
2. This unit has an optional AC Power adaptor which is available as an accessory. Only use AC adaptor with below specification (not included).  
Input: 100-240VAC 50/60Hz 0.2A Max  
Output: 6V  $\equiv$  1000mA  
(Conforms to UL certificate)  
Right picture is the hole in for power adaptor.

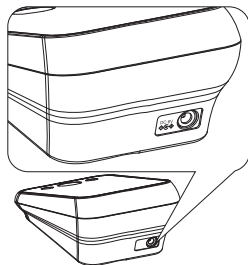


### ⚠ CAUTION


In order to get the best effect and protect your monitor, please use the right battery and special power adaptor.

## ♥ Installing and Replacing the Batteries

1. Slide off the battery cover.
2. Install the batteries by matching the correct polarity, as shown.
3. Replace the cover.

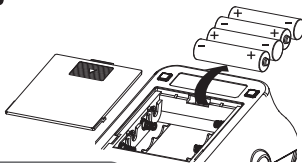

Replace the batteries whenever the below happen

- The  shows
- The display dims
- The display does not light up


### ⚠ CAUTION

- Remove batteries if the device is not likely to be used for some time.
- The old batteries are harmful to the environment, so please do not dispose with other daily trash.
- Remove the old batteries from the device and follow your local recycling guidelines.
- Do not dispose of batteries in fire. Batteries may explode or leak.


## ♥ Setting Date, Time and Measurement Unit

It is important to set the clock before using your blood pressure monitor, so that a time stamp can be assigned to each record that is stored in the memory. (year :2018—2058 time format:12 H)


1. When the unit is off, hold "SET" for 3 seconds to enter the mode for year setting.

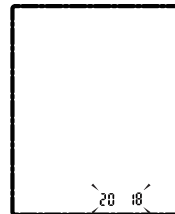


2. Press the "MEM" to change the [YEAR]. Each press will increase the numeral by one in a cycling manner.

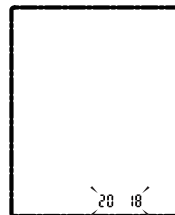


3. When you get the right year, press "SET" to confirm, and it will divert to the [MONTH] setting.

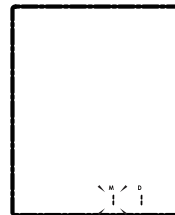


4.Repeat steps 2 and 3 to set the [MONTH] and [DAY].

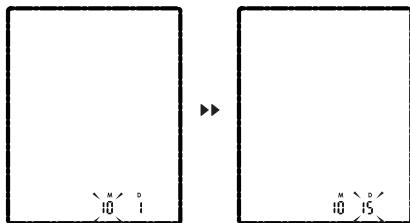

5.Repeat steps 2 and 3 to set the [HOUR] and [MINUTE].

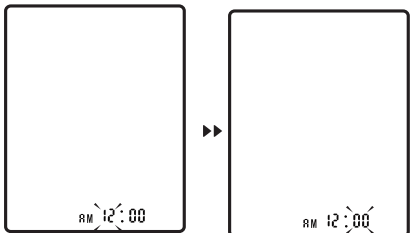

6.Repeat steps 2 and 3 to set the [MEASUREMENT UNIT].

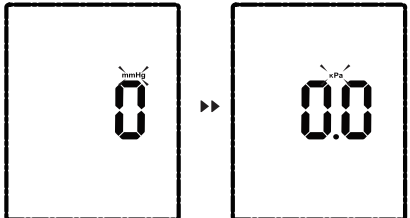

7.After the [MEASUREMENT UNIT] is set,the LCD will display “dOnE”, and then turn off.



## ♥ Tie the Cuff

1. Remove all jewelry, such as watches and bracelets from your left arm.  
Note: If your doctor has diagnosed you with poor circulation in your left arm, use your right arm.

2. Roll or push up your sleeve to expose the skin. Make sure your sleeve is not too tight.



3. Hold your arm with your palm facing up and tie the cuff on your upper arm, then position the tube off-center toward the inner side of arm in line with the little finger. Or position the artery mark  $\Phi$  over the main artery (on the inside of your arm). Note: Locate the main artery by pressing with 2 fingers approximately 2 cm above the bend of your elbow on the inside of your left arm. Identify where the pulse can be felt the strongest. This is your main artery.

4. The cuff should be snug but not too tight. You should be able to insert one finger between the cuff and your arm.

5. Sit comfortably with your tested arm resting on a flat surface. Place your elbow on a table so that the cuff is at the same level as your heart. Turn your palm upwards. Sit upright in a chair, and take 5-6 deep breaths.

6. Helpful tips for Patients, especially for Patients with Hypertension:

- Rest for 5 minutes before first measurement.
- Wait at least 3 minutes between measurements. This allows your blood circulation to recover.
- Take the measurement in a silent room.
- The patient must relax as much as possible and do not move and talk during the measurement procedure.
- The cuff should maintain at the same level as the right atrium of the heart.
- Please sit comfortably. Do not cross your legs and keep your feet flat on the ground.
- Keep your back against the backrest of the chair.
- For a meaningful comparison, try to measure under similar conditions. For example, take daily measurements at approximately the same time, on the same arm, or as directed by a physician.

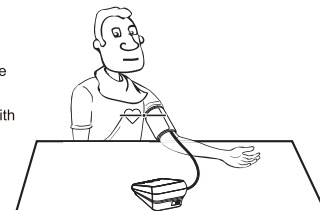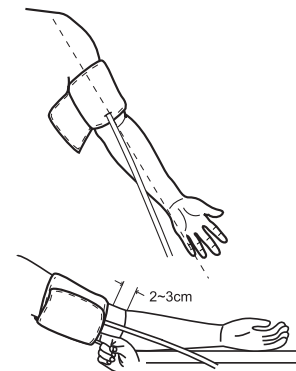


## ♥ Start the Measurement

Before you start the measurement, please press the SET button to choose either User 1 or User 2 as the User ID. When the desired User ID is shown, press START/STOP button to confirm the User ID.

1. After selecting the user, press the "START/STOP" to start measurement, and it will finish the whole measurement for the selected user.

Take User 1 for example:

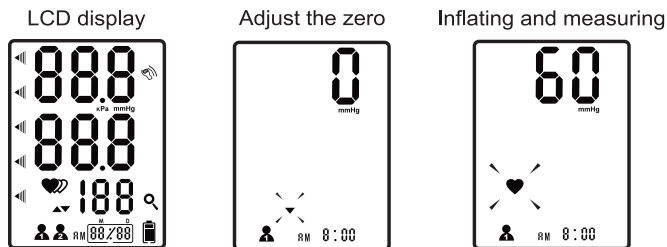

Display and save the results. The corresponding backlight shows according to the grade of blood pressure.




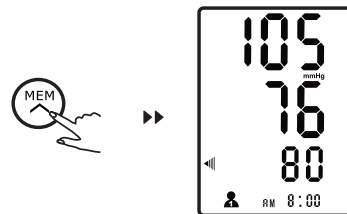
2. Press "START/STOP" button to power off, otherwise it will turn off within 1 minute.

### Tips:

- A. You can press "START/STOP" button at any time to stop measuring.
- B. Maximum 60 records are recorded for both for USER 1 and USER 2.
- C. If the measurement result is out of the measurement range (SYS: 60mmHg -230mmHg; or DIA:40mmHg - 130mmHg; or Pulse: 40-199 pulse/minute), the LCD will display "out". Consult the Error code on page 18 for more information and further instructions.


## ♥ Recall the Records

1. When the monitor is off, please press the "MEM" to show the average value of the latest three records. If there are less than 3 record entries, it will display the latest record instead and the icon  will not appear. Device displays different colors at different levels when the monitor is under the memory mode. Consult the standard blood pressure classification on pages 16 for more information and further instructions.



2. Press MEM button or SET button to rotate the records. Up to 60 records will be stored under each user ID.

Date,time will be shown under the Pulse by turns.

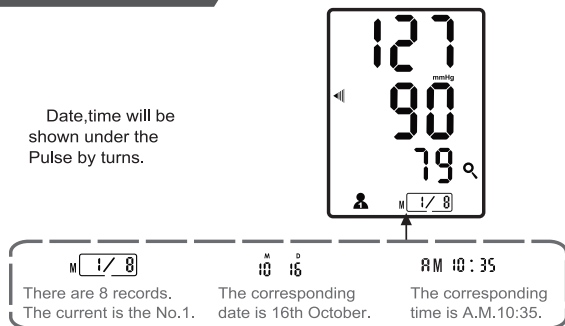

3. If you want to check another user's records, press STAR/SOP button to turn off the monitor when it is in the memory recall mode. Then press SET button, the user icon will be shown, press "SET" button to select the desired user ID, press "MEM" button to review the selected user's records.



4. Press the START/STOP button to turn off the monitor. Otherwise, the monitor will shut off within 1 minute after last operation.


### ⚠ CAUTION

The most recent record (1) is shown first. Each new measurement is assigned to the first (1) record. All other records are pushed back one digit (e.g., 2 becomes 3, and so on), and the last record (60) is dropped from the list.

## ♥ Delete the Records

If you did not get the correct measurement, you can delete all results or the latest result for the selected user by the following steps below.

1. If you want to delete all the results, hold "MEM" for 3 seconds when the monitor is in the memory mode, then "dEL ALL" will show.

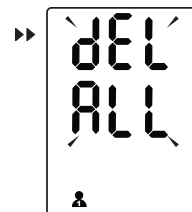


OR


- If you want to delete the latest result, hold "SET" for 3 seconds when the monitor is in the memory mode, then "dEL OnE" will show.



2. Press "MEM" to confirm deleting and the monitor will turn off.

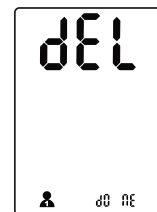

3. If you don't want to delete the records, press "STRAT/STOP" button or "SET" button to escape.

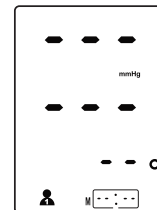

4. If there is no record, the right display will show when recalling the record.



♥ Tips for Measurement

Measurements may be inaccurate if taken in the following circumstances.



wait at least 1 hour after dinner or drinking



Wait at least 20 minutes after taking a bath



In a very cold environment



Immediate measurement after tea, coffee, smoking



When talking or moving your fingers




When you want to discharge urine

♥ Maintenance

In order to get the best performance, please follow the instructions below.



Put in a dry place and avoid the sunshine

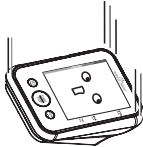

Avoid intense shaking and collisions



Use wet cloths to remove dirt



Avoid touching water, clean it with a dry cloth in case.




Avoid dusty and unstable temperature environment



Avoid washing the cuff

♥ What are systolic pressure and diastolic pressure?

When ventricles contract and pump blood out of the heart, the blood pressure reaches its maximum value in the cycle, which is called systolic pressure. When the ventricles relax, the blood pressure reaches its minimum value in the cycle, which is called diastolic pressure.

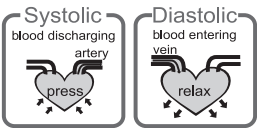

♥ What is the standard blood pressure classification?

The chart on the right is the standard blood pressure classification published by American Heart Association (AHA).

AHA Home Guideline for Upper Limit of Normal BP

SYS	135 mmHg
DIA	85 mmHg

This chart reflects blood pressure categories defined by American Heart Association.					
Backlight on BPM	AHA Color guide line	Blood Pressure Category	Systolic mmHg (upper#)		Diastolic mmHg (lower#)
Blue	Green	Normal	less than 120	and	less than 80
Blue	Green	Prehypertension	120-129	and	less than 80
Orange	Yellow	High Blood Pressure (Hypertension) Stage 1	130-139	or	80-89
Red	Orange	High Blood Pressure (Hypertension) Stage 2	140 or higher	or	90 or higher
Red	Red	Hypertensive Crisis (Emergency care needed)	Higher than 180	and/or	Higher than 120

⚠ CAUTION

Only a physician can tell your normal BP range. Please contact a physician if your measuring result falls out of the range. Please note that only a physician can tell whether your blood pressure value has reached a dangerous point.

♥ Irregular Heartbeat Detector


An irregular heartbeat is detected when a heartbeat rhythm varies while the device is measuring systolic pressure and diastolic pressure. During each measurement, blood pressure monitor will keep a record of all the pulse intervals and calculate the average value of them. If there are two or more pulse intervals, the difference between each interval and the average is more than the average value of  $\pm 25\%$ , or there are four or more pulse intervals, the difference between each interval and the average is more than the average value of  $\pm 15\%$ , then the irregular heartbeat symbol will appear on the display with the measurement result to indicate this abnormality.

⚠ CAUTION

The appearance of the IHB icon indicates that a pulse irregularity consistent with an irregular heartbeat was detected during measurement. Usually this is NOT a cause for concern. However, if the symbol appears often, we recommend you seek medical advice. Please note that the device does not replace a cardiac examination, but serves to detect pulse irregularities at an early stage.

♥ Why does my blood pressure fluctuate throughout the day?

1. Individual blood pressure varies multiple times everyday. It is also affected by the way you tie your cuff and your measurement position, so please take the measurement under the same conditions.
- 2.If the person takes medicine, the pressure will vary more.
- 3.Wait at least 3 minutes for another measurement.

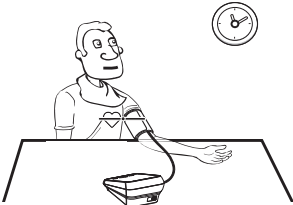

♥ Why do I get a different blood pressure at home compared to the hospital?


The blood pressure is different even throughout the day due to weather, emotion, exercise etc. Also, there is the “white coat” effect, which means blood pressure usually increases in clinical settings.

What you need to pay attention to when you measure your blood pressure at home:


- If the cuff is tied properly.
- If the cuff is too tight or too loose.
- If the cuff is tied on the upper arm.
- If you feel anxious.
- Taking 2-3 deep breaths before beginning will be better for measuring.
- Advice: Relax yourself for 4-5 minutes until you calm down.

♥ Is the result the same if measuring on the right arm?

It is ok for both arms, but there will be some different results for different people. We suggest you measure the same arm every time.

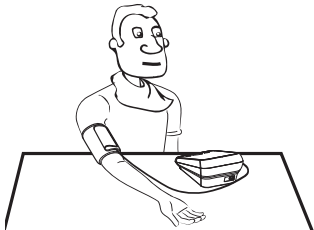


This section includes a list of error messages and frequently asked questions for problems you may encounter with your blood pressure monitor. If the products not operating as you think it should, check here before arranging for servicing.

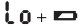PROBLEM	SYMPTOM	CHECK THIS	REMEDY
No power	Display is dim or will not light up.	Batteries are exhausted.	Replace with new batteries
		Batteries are inserted incorrectly.	Insert the batteries correctly
		AC adaptor is inserted incorrectly.	Insert the AC adaptor tightly
Low batteries	 Show on the display	Batteries are low.	Replace with new batteries
Error message	E 3 shows	The cuff is not secure.	Refasten the cuff and then measure again.
	E 10 or E11 shows	The monitor detected motion while measuring.	Movement can affect the measurement. Relax for a moment and then measure again.
	E20 shows	The measurement process does not detect the pulse signal.	Loosen the clothing on the arm and then measure again
	E21 shows	The treatment of the measurement failed.	Relax for a moment and then measure again.
	EExx, shows on the display.	A calibration error occurred.	Retake the measurement. If the problem persists, contact the retailer or our customer service department for further assistance. Refer to the warranty for contact information and return instructions.
Warning message	"out " shows	Out of measurement range	Relax for a moment. Refasten the cuff and then measure again. If the problem persists, contact your physician.

Power supply	Battery powered mode: 6VDC 4*AA batteries AC adaptor powered mode: (INPUT: 100-240VAC 50/60Hz 0.2A Max OUTPUT: 6V--- 1000mA)(Not Included)
Display mode	Digital LCD V.A.60mm*80mm (2.36"*3.15")
Measurement mode	Oscillographic testing mode
Measurement range	Rated cuff pressure: 0mmHg~299mmHg(0kPa ~ 39.9kPa) Measurement pressure: SYS: 60mmHg~230mmHg (8.0kPa~30.7kPa) DIA: 40mmHg~130mmHg (5.3kPa~17.3kPa) Pulse value: (40-199)beat/minute
Accuracy	Pressure: 5°C~40°C(41°F~104°F) within ±3mmHg(0.4kPa) pulse value: ±5%
Normal working condition	A temperature range of :+5°C to +40°C A relative humidity range of 15% to 90%, non-condensing, but not requiring a water vapour partial pressure greater than 50 hPa An atmospheric pressure range of : 700 hPa to 1060 hPa
Storage & transportation condition	Temperature:-20°C to +60°C A relative humidity range of ≤ 93%, non-condensing, at a water vapour pressure up to 50hPa
Measurement perimeter of the upper arm	About 22cm~32cm ( 8 1/4"~12 1/2" ) or 22cm~42cm (8 1/4"~16 1/2" )
Net Weight	Approx.262g( 9.24oz)(Excluding the batteries)
External dimensions	Approx.102mm*143mm*73mm(4.02"*5.63"*2.87")
Attachment	4*AA batteries, one storage bag, user manual
Mode of operation	Continuous operation
Degree of protection	Type BF applied part
Protection against ingress of water	IP21: This device is protected against solid foreign objects of 12.5mm and greater, and protected against vertically falling water drops.
Software Version	A01

## ♥ Authorized Components

1. please use the  authorized adapter. (Not Included)

- 2.Storage bag.



### Adapter

Input: 100-240VAC 50/60Hz 0.2A Max  
Output: 6V  $\equiv$  1000mA  
( BLJ06L060100P-U )

## ♥ Contact Information

For more information about our products, please visit [JHealthcare-Manager.com](http://JHealthcare-Manager.com), or call toll-free at 1-855-822-6999 M-F 9 a.m.-5 p.m. CST. We attend to all questions, concerns, and guidance about our resources.

## ♥ Contraindications

- 1.The device should not be used by any person who can possible be or is pregnant.
- 2.The device is not suitable for use on patients with implanted, electrical devices, such as cardiac pacemakers, defibrillators.

## ♥ Measurement Principle

This product uses the Oscillometric Measuring Method to detect blood pressure. Before every measurement, the unit establishes a "zero point" equivalent to the atmospheric pressure. Then it starts inflating the cuff. Meanwhile, the unit detects pressure oscillation generated by beat-to-beat pulsation, which is used to determine the systolic pressure and diastolic pressure as well as pulse rate. Systolic pressure is the maximum pressure your heart exerts while beating (the top number), and diastolic pressure is the amount of pressure in your arteries between beats (the bottom number). Pulse rate (also known as heart rate) is the number of times your heart beats per minute.

## ♥ Complied Standards List

<b>Risk management</b>	EN ISO 14971:2012 / ISO 14971:2007 Medical devices - Application of risk management to medical devices
<b>Labeling</b>	EN ISO 15223-1:2016 / ISO 15223-1:2016 Medical devices. Symbols to be used with medical device labels, labelling and information to be supplied. Part 1 : General requirements
<b>User manual</b>	EN 1041:2008 +A1:2013 Information supplied by the manufacturer of medical devices
<b>General Requirements for Safety</b>	EN 60601-1:2006+A1:2013/ IEC 60601-1:2005+A1:2012 Medical electrical equipment - Part 1: General requirements for basic safety and essential performance EN 60601-1-11:2015/ IEC 60601-1-11:2015 Medical electrical equipment - Part 1-11: General requirements for basic safety and essential performance - Collateral standard: Requirements for medical electrical equipment and medical electrical systems used in the home healthcare environment
<b>Electromagnetic compatibility</b>	EN 60601-1-2:2015/ IEC 60601-1-2:2014 Medical electrical equipment - Part 1-2: General requirements for basic safety and essential performance - Collateral standard: Electromagnetic disturbances - Requirements and tests
<b>Performance requirements</b>	EN ISO 81060-1:2012 Non-invasive sphygmomanometers - Part 1: Requirements and test methods for non-automated measurement type EN 1060-3:1997+A2:2008 Non-invasive sphygmomanometers - Part 3: Supplementary requirements for electro-mechanical blood pressure measuring systems IEC 80601-2-30:2009+A1:2013 Medical electrical equipment- Part 2-30: Particular requirements for the basic safety and essential performance of automated non-invasive sphygmomanometers
<b>Clinical investigation</b>	EN 1060-4:2004 Non-invasive sphygmomanometers - Part 4: Test procedures to determine the overall system accuracy of automated non-invasive sphygmomanometers ISO 81060-2:2013 Non-invasive sphygmomanometers - Part 2: Clinical validation of automated measurement type
<b>Usability</b>	EN 60601-1-6:2010+A1:2015/IEC 60601-1-6:2010+A1:2013 Medical electrical equipment - Part 1-6: General requirements for basic safety and essential performance - Collateral standard: Usability IEC 62366-1:2015 Medical devices - Part 1: Application of usability engineering to medical devices
<b>Software life-cycle processes</b>	EN 62304:2006/AC: 2008 / IEC 62304: 2006+A1:2015 Medical device software - Software life-cycle processes
<b>Bio-compatibility</b>	ISO 10993-1:2009 Biological evaluation of medical devices- Part 1: Evaluation and testing within a risk management process ISO 10993-5:2009 Biological evaluation of medical devices - Part 5: Tests for in vitro cytotoxicity ISO 10993-10:2010 Biological evaluation of medical devices - Part 10: Tests for irritation and skin sensitization

♥ FCC Statement

This device complies with Part 15 of the FCC Rules. Operation is subject to the following two conditions: (1) this device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.

♥ EMC Guidance

The ME EQUIPMENT or ME SYSTEM is suitable for home healthcare environments

Warning: Don't near active HF surgical equipment and the RF shielded room of an ME system for magnetic resonance imaging, where the intensity of EM disturbances is high.

Warning: Use of this equipment adjacent to or stacked with other equipment should be avoided because it could result in improper operation. If such use is necessary, this equipment and the other equipment should be observed to verify that they are operating normally.

Warning: Use of accessories, transducers and cables other than those specified or provided by the manufacturer of this equipment could result in increased electromagnetic emissions or decreased electromagnetic immunity of this equipment and result in improper operation."

Warning: Portable RF communications equipment (including peripherals such as antenna cables and external antennas) should be used no closer than 30 cm (12 inches) to any part of the equipment EBP-095L, including cables specified by the manufacturer. Otherwise, degradation of the performance of this equipment could result.

Technical description:

- 1, all necessary instructions for maintaining BASIC SAFETY and ESSENTIAL PERFORMANCE with regard to electromagnetic disturbances for the excepted service life.
- 2, Guidance and manufacturer's declaration -electromagnetic emissions and Immunity

Table 1

Guidance and manufacturer's declaration - electromagnetic emissions	
Emissions test	Compliance
RF emissions CISPR 11	Group 1
RF emissions CISPR 11	Class [ B ]
Harmonic emissions IEC 61000-3-2	Class A
Voltage fluctuations/ flicker emissions IEC 61000-3-3	Comply

Table 2

Guidance and manufacturer's declaration – electromagnetic immunity		
Immunity Test	IEC 60601-1-2 Test level	Compliance level
Electrostatic discharge (ESD) IEC 61000-4-2	±8 kV contact ±2 kV, ±4kV, ±8 kV, ±15 kV air	±8 kV contact ±2 kV, ±4kV, ±8 kV, ±15 kV air
Electrical fast transient/burst IEC 61000-4-4	±2 kV for power supply lines ±1 kV signal input/output 100 kHz repetition frequency	±2 kV for power supply lines ±1 kV signal input/output 100 kHz repetition frequency
Surge IEC61000-4-5	±0.5 kV, ±1 kV differential mode ±0.5 kV, ±1 kV,±2 kV common mode	±0.5 kV, ±1 kV differential mode ±0.5 kV, ±1 kV,±2 kV common mode

Voltage dips, short interruptions and voltage variations on power supply input lines IEC 61000-4-11	0 % U <sub>T</sub> ; 0,5 cycle. At 0°, 45°, 90°, 135°, 180°, 225°, 270° and 315°. 0 % U <sub>T</sub> ; 1 cycle and 70 % U <sub>T</sub> ; 25/30 cycles; Single phase: at 0°. 0 % U <sub>T</sub> ; 250/300 cycle	0 % U <sub>T</sub> ; 0,5 cycle. At 0°, 45°, 90°, 135°, 180°, 225°, 270° and 315°. 0 % U <sub>T</sub> ; 1 cycle and 70 % U <sub>T</sub> ; 25/30 cycles; Single phase: at 0°. 0 % U <sub>T</sub> ; 250/300 cycle
Power frequency magnetic field IEC 61000-4-8	30 A/m 50Hz/60Hz	30 A/m 50Hz/60Hz
Conducted RF IEC61000-4-6	3 V 0,15 MHz – 80 MHz 6 V in ISM and amateur radio bands between 0,15 MHz and 80 MHz 80 % AM at 1 kHz	3 V 0,15 MHz – 80 MHz 6 V in ISM and amateur radio bands between 0,15 MHz and 80 MHz 80 % AM at 1 kHz
Radiated RF IEC61000-4-3	10 V/m 80 MHz – 2,7 GHz 80 % AM at 1 kHz	10 V/m 80 MHz – 2,7 GHz 80 % AM at 1 kHz
NOTE U <sub>T</sub> is the a.c. mains voltage prior to application of the test level.		

810	800-960	GSM 800/900, TETRA 800, iDEN 820, CDMA 850, LTE Band 5	Pulse modulation b) 18Hz	2	0.3	28
870						
930						
1720	1700-1990	GSM 1800; CDMA 1900; GSM 1900; DECT; LTE Band 1, 3, 4, 25; UMTS	Pulse modulation b) 217Hz	2	0.3	28
1845						
1970						
2450	2400-2570	Bluetooth, WLAN, 802.11 b/g/n, RFID 2450, LTE Band 7	Pulse modulation 217 Hz	2	0.3	28
5240	5100-5800	WLAN 802.11 a/n	Pulse modulation 217 Hz	0.2	0.3	9
5500						
5785						

Table 3

Guidance and manufacturer's declaration - electromagnetic Immunity							
Radiated RF IEC61000-4-3 (Test specifications for ENCLOSURE PORT IMMUNITY to RF wireless communications equipment)	Test Frequency (MHz)	Band (MHz)	Service	Modulation	Modulation (W)	Distance (m)	IMMUNITY TEST LEVEL (V/m)
	385	380-390	TETRA 400	Pulse modulation b) 18Hz	1.8	0.3	27
	450	430-470	GMRS 460, FRS 460	FM c) ± 5kHz deviation 1kHz sine	2	0.3	28
	710	704-787	LTE Band 13, 17	Pulse modulation b) 217Hz	0.2	0.3	9
	745						
	780						