


Thank you for purchasing Easy@Home Water Flosser.

Water Flosser provides a unique combination of water pressure and pulsation to clean deep between teeth and the gum, removing debris of food and harmful bacteria that conventional brushing or flossing can't reach. It also massages and stimulates the gum to improve blood circulation and helps in keeping the gums strong and healthy. It is also ideal for cleaning around dental bridges, braces and crowns. Enjoy oral SPA at home with Easy@Home Water Flosser!

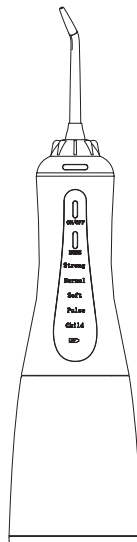

## Water Flosser

User Manual

**EWF003**

Thank you for choosing Easy@Home Water Flosser. Make your teeth healthier, and your smile brighter!

Safety Warning .....	01
Precautions .....	02
Product Description .....	03
Operating Procedures .....	04
Usage Instructions .....	06
Recommended Techniques .....	07
Storage Guidelines .....	08
Battery Maintenance And Handling .....	09
All The Accessories .....	10
Common Problems And Solutions .....	11
Contents .....	11
Expiry Date / Service Life / Warranty .....	11
Storage And Transportation Condition .....	12
Bibliography / Complied Standards .....	12
FCC And EMC Information .....	13
Explanation Of Symbols .....	14

## Contraindications

**Read this instruction book thoroughly before use.**


- ◆ Do not immerse the flosser in water or any other liquid.
- ◆ Do not use it while bathing or showering.
- ◆ Consult your dentist any time if you experience discomfort while using this water flosser.

## WARNING

**This product can be used by children ages 8 and above and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge only if they are under supervision or have been given instructions on the proper use of the product in a safe way and understand the hazards involved.**


- ◆ Water Flosser shall be used only with cold or warm water, or specific solutions as defined in the instruction manual.
- ◆ Keep away from small children.
- ◆ This product contains built-in batteries that are non-replaceable.

- ◆ If the product is opened/disassembled for ANY reason, warranty is VOID.
- ◆ Do not use the product if you see a damaged charging cord or defect at any part of the product.
- ◆ Do not direct water under the tongue, into the ear, nose, or other delicate areas. This product is capable of producing pressures that may cause serious damage in these areas.
- ◆ Do not direct water into the nose or ear. The potentially deadly amoeba, *Naegleria fowleri*, may be present in some tap water or unchlorinated well water and may be fatal if directed into these areas.
- ◆ Use this product only as indicated in these instructions or as recommended by your dental professional.
- ◆ Magnetic field: Pacemakers and similar implanted medical devices may be affected by the magnet produced by this product while running. If you have an implanted medical device, please consult your physician before purchasing and using.
- ◆ Do not drop or insert any foreign object into any opening of the product.
- ◆ Fully charge the device before using.
- ◆ Do not use this product outdoors or operate where aerosol (spray) products are being used or where oxygen is being administered.
- ◆ Children should be supervised to ensure that they do not play with this product.

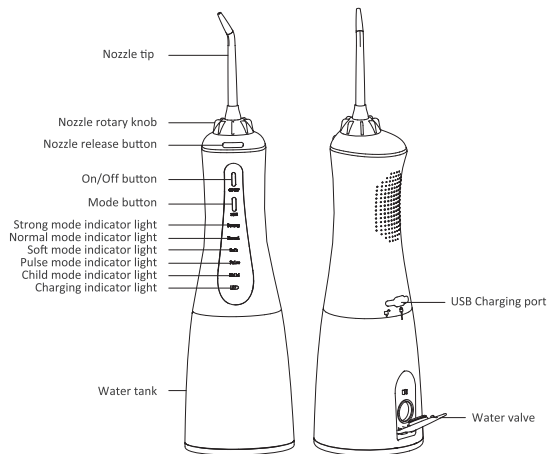


## Preparation Work.

Important! Before using for the first time, please fully charge the device. During charging, the product CAN NOT be used.

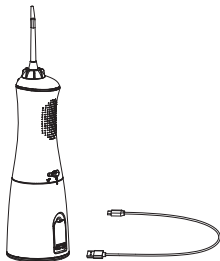


## Insert and remove tips


Insert the tip into the center of the knob until there is "click" sound.

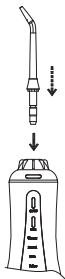

Open the water tank lid and pour water into the water tank;  
Note: please ensure that when pouring water into the water tank, the device is switched off.

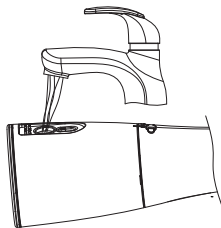

To remove the tip, press the nozzle release button.

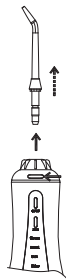


**This product has a memory function. The mode remains the same as the last mode used when switched on again.**

**Please select mode before switching the product on.**


- ◆ The charging indicator light stays red when charging, and stays green when it is fully charged.
- ◆ When charging indicator light becomes red and flashes, it means the battery is low and recharging is required.
- ◆ This product can NOT be used during charging.
- ◆ With the ergonomically designed and 360° rotated tip, you can easily control water flow to reach all areas of the mouth.



Direct the jet stream at a 90° angle to your gumline.

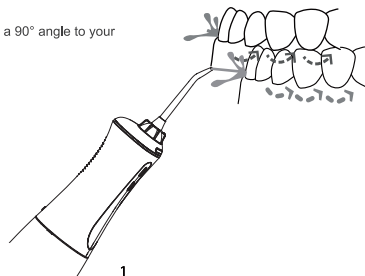

Close lips to avoid splashing but **allow** water to flow freely from mouth into the sink. For best result, start in the molar area (back teeth) working toward the front teeth. Glide jet tip along gumline and pause briefly between teeth. Continue until **all** areas around and between the teeth have been cleaned.

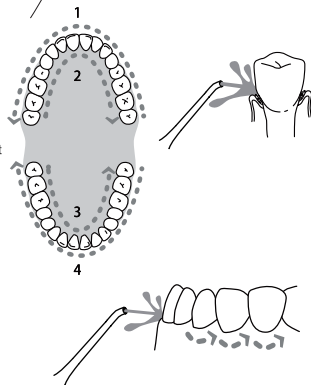


1

Please switch off the power, as well as empty the water tank.




2

Use both hands rotating downwards in clockwise direction as indicated in Diagram 2 and 3 to remove the water tank from the body.

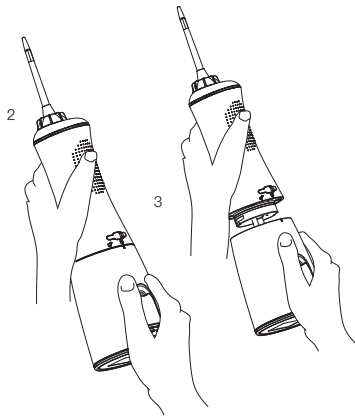

3

Removed nozzle tips can be stored in the water tank, as illustrated in Diagram 4.

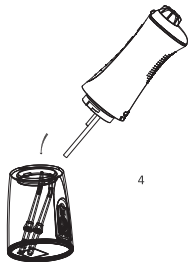






## Battery Maintenance And Handling

### Battery Maintenance

- ◆ For the first use, the product needs to be fully charged, for approx. 5 hours.
- ◆ If the product is to be stored for a long time, please fully charge it in advance.
- ◆ In order to extend battery lifespan, recharge it in time when low power indication is seen.

### Battery Handling

The product has built-in lithium battery which cannot be replaced. To protect the environment, do not throw away the product with other normal household waste at the end of its life. Deliver it to an official collection point for recycling.

Pictures	Description	Function
	Jet Tips	For general use. Cleans deep between teeth and below the gumline. Clinically proven to be twice as effective as string floss for improving gum health. Replace every 3 months
	Orthodontic Tip	For orthodontic appliances and dental braces. Clinically proven to remove three times as much plaque as traditional string floss and significantly improve gum health for orthodontic patients. The innovative tapered brush removes hard-to-reach plaque around orthodontic brackets and wires. It also helps flush out bacteria and food debris from around the teeth and under the gums. Replace every 3 months
	Pocket Tip	For periodontal pockets and furcations. A soft rubber tip designed for gentle, low-pressure delivery of therapeutic rinses below the gumline. Clinically proven to access 90% of a 6mm pocket. Replace every 3 months
	Plaque tip	The Tips use for dental areas that are susceptible to plaque accumulation, including implants, crowns, bridges, and permanent retainers. Replace every 3 months
	Tongue Cleaner	For fresh breath. With a gentle scraping and flushing action, the Tongue Cleaner removes bacteria and sulfur compounds that reside on the tongue and cause bad breath. Daily use can help prevent halitosis. Replace every 6 months

Note:  
All the accessories we quote

Problems	Possible Causes	Solutions
Leakage of water at connection parts	1. The jet nozzle or the water tank are not set in position. 2. The o-ring seal at the connection part is missing or damaged.	1. Set them in position. 2. Contact the seller for assistance.
Inadequate water pressure	Low power	Charge immediately
Unable to switch on	Low power	Charge immediately
Motor is running but water is not sprayed out	The nozzle is clogged or the water valve seal is not tight	Dredge the nozzle, make the water valve seal tight

## Contents

1 x Water Flosser	1 x Tongue Cleaner
2 x Jet Tips	1 x Storage bag
1 x Orthodontic Tip	1 x USB charging cable
1 x Pocket Tip	1 x User manual
1 x Plaque tip	

## Expiry date / Service life / Warranty

Expiry date: 3 years

Warranty: 1 year

**Working environment**

Temperature: 5~40°C

Humidity: 15% to 85%RH

Atmospheric pressure: 86Kpa to 106 Kpa

**Transport & storage environment**

Temperature: -20~55°C

Humidity: ≤85%RH

Atmospheric pressure: 86Kpa to 106 Kpa

**STANDARDS USED(DATE) :**

EN 55014-1:2017 Electromagnetic compatibility - Requirements for household appliances, electric tools and similar apparatus Part 1: Emission (CISPR 14-1:2016 + COR1:2016)

EN 55014-2:2015 EN 55014-2:2015 Electromagnetic compatibility - Requirements for household appliances, electric tools and similar apparatus - Part 2: Immunity - Product family standard

IEC 60335-2-52 Household and similar electrical appliances – Safety – Part 2-52: Particular requirements for oral hygiene appliances

EMC Emission			
Standard	Test Item	Limit	Result
FCC Part 15 Subpart B ANSI C63.4: 2014	Conducted Emission	Class B	N/A
	Radiated Emission	Class B	PASS

Note: 1. "N/A" denotes test is not applicable in this test report







2. For client's request and manual description, the test will not be executed.

3. This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to part 15 of the FCC Rules.

Standard	Test Item	Limit	Judgment
BS EN 55014 -1:2017/A11:2020	Conducted Emission	Class B	N/A
	Disturbance Power(30 ~300MHz)	—	N/A
	Radiated Emission 30MHz to 1000MHz	Class B	PASS

EMC Immunity			
Section BS EN 55014 -2:2015	Test Item	Performance Criteria	Judgment
BS EN 61000 -4-2:2009	Electrostatic Discharge	B	PASS
BS EN 61000 -4-3:2006 +A1:2008+A2:2010	RF Electromagnetic Field	A	PASS



	WEEE (Waste Electrical and Electronic Equipment)		Recyclable
	CE marking of conformity		Federal Communications Commission
	Manufacturing date		Batch code

**Healthcare-Manager.com**

Questions or comments?

Please call toll-free:

1-855 822 -6999 M-F 9 a.m-5 p.m. CST

E-MAIL: [service@healthcare-manager.com](mailto:service@healthcare-manager.com)

**Manufactured for**

**Easy Healthcare Corporation**

360 Shore Dr. Unit B, Burr Ridge, IL USA 60527

Made in China

V01: 20220407