


HEATED EYE MASK




Note: The eye mask will emit heat 3-5 minutes after turned on

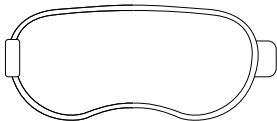


Made in China



Using Instruction

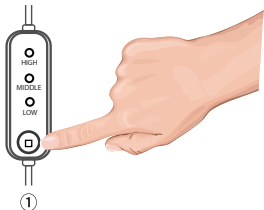

Power / Temperature Button

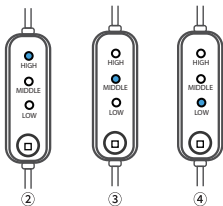

② HIGH: 55 °C

③ MIDDLE: 45 °C

④ LOW: 40 °C

- 1.Connect the product to the power supply using the USB cable.
 - 2.Press the button to turn on the product, and the default level is “HIGH”.
 - 3.Press the button for the second time to activate “MIDDLE” level.
 - 4.Press the button for the third time to activate “LOW” level.
 - 5.When the product works at “LOW” level, press the button once to turn off it.
- *The product will shut off automatically after 30 minutes of continuous use.

NOTE:

1. The eye mask can be washed, but it's not machine washable. DO NOT wash violently.
2. Do not use in a humid environment such as a bathroom.
3. Do not use if you have eye surgery or bleeding tendency and other blood disorders.
4. DO NOT use in the case of osteoporosis or soft tissue disorder.
5. Please stop using if there is any discomfort during use.
6. Children under the age of 6 are forbidden to use and children aged 6-12 should use under the guidance and care.
7. Do not use sharp objects to scratch the product.

Welcome to MOUNTRAX

Want to explore more with MOUNTRAX? You may visit our website.

About customer service, please contact our customer support email.

Welcome to share your experience with MOUNTRAX, leave your story on our Facebook!



www.mountrax.net



support@mountrax.net



facebook.com/MountraxPage

SPECIFICATIONS:

Model: EM-004

Rated Voltage: 5V

Current: 1A-2A (Max.)

Max. Temperature: 55°C

Timer: 30 minutes

Package Contents: USB cable x 1, Heated eye mask x 1