

biroller® **G5**



User Manual

Please read the manual before use

Product description	2
Parts name	3
Operation	4
Operation	5
Operation	6
How to charge the battery	7
Specification	8
Warranty	9
Precaution	10

The device consists of an electric handpiece and a microcurrent roller.

This manual will explain how to use each component in detail.


The device combines the technology of microneedling with microcurrent.

Advantage of microneedling:


Improve the appearance of the skin and reduce signs of aging. It is intended to stimulate the skin with micro-needles for collagen and elastin production. The micro-channels created on surficial skin by the device also help with the absorption of serums.

Advantage of microcurrent:

Enhance the activity of cells, promote muscle movement and restore skin elasticity. Accelerate the blood circulation of capillaries

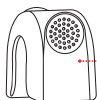

Protective cap

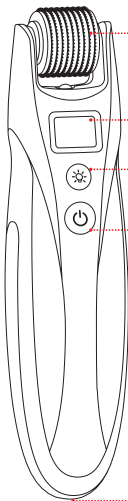

Microcurrent Roller


Display Screen

Light button
Vibration button

ON / OFF
Intensity button


Power Interface

1. Take out needle cartridge from the sterilized packaging.
〔Fig1〕

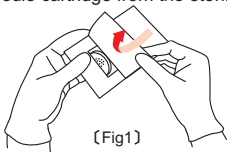


2. Please hold the side with "SKIN" mark upward.
And then assemble the roller with protective cap to the connector of handpiece. 〔Fig2〕

Attention: Please sterilize the roller with 75% alcohol for few minutes before assemble the roller to the handpiece.

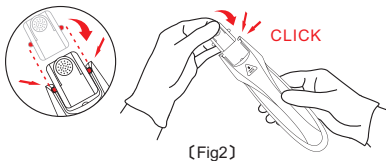

3. Remove the protective cap from the roller.
It is ready to use now. 〔Fig3〕

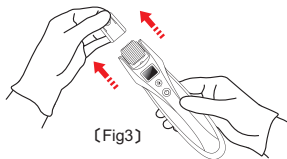

4. Press and hold the power button for 3 seconds to start or shut down the device.
5. 5 different levels for the intensity of microcurrent. Use suitable level for intensity according to the treatment and skin type. Change the intensity of microcurrent by pressing the intensity button. The intensity of microcurrent varies with each press. [Fig4]



[Fig4]

6. Press the light button to start or shut down the light and vibration. The light color and vibration changes each time the button is pressed [Fig5]



[Fig5]



Blue light (Weak vibration)


It is good for sterilization and improve the acne.



Red light (Strong vibration)

Reduce inflammation and stimulate the production of collagen.

7. Hold the side with "SKIN" mark toward the skin and then roll on the skin for treatment. [Fig6]

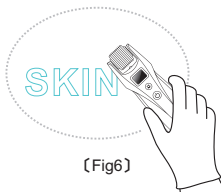


[Fig6]



Ensure this side facing Skin while working

How to remove the microcurrent roller

1. Assemble the protective cap back to the roller [Fig7]

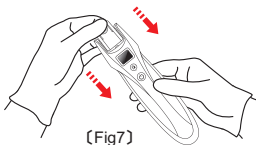

[Fig7]

2. Disassemble the roller with protective cap from the handpiece [Fig8]



[Fig8]

Wireless Model

This device is powered by a built-in battery. It can work wirelessly. The power cable and the adapter are used to charge the device when the battery is low.


Connect the power cable with adapter to the power interface of device and plug the adapter into a live AC outlet to charge the device. Charging will complete in about 2 hours.


Charging Status

The number of speed flashes slowly when the battery's capacity 10% left .

The pen will turn off automatically after the number of speed flashes quickly for 1 minute.

The pen shows “  ” when restarting after the low battery automatic shutdown.

The pen show spinning “  ” while charging .

The pen show “  ” when charging is completed.

Model name: bio roller G5

Input: 100-240V 50/60Hz

Output: 5V \equiv 2A

Package Size: 123*168*56 mm

Package Weight: 350g

Packing List

Item	Quantity
Handpiece	1 pc
Massage roller	2 pcs
Manual	1 pc
Power cable	1 pc
Adapter	1 pc

The warranty period of the device is 12 months. We offer repair service from the purchase date. If any of our products has damages or natural trouble, we provide replacements or repair service.

General maintenance

Make sure to inspect the device and adapter before use.

Ensure the device is turned off when not in use. Wipe its surface with a disinfecting wipe.

Pay attention to clean areas such as crevices, seams, and areas where the microcurrent roller attaches to the device.

Keep it away from direct sunlight or high temperature or humidity.

Do not leave the measurement, repair or alter of the device to non-professional staff.

When you use the device which has not been used for a long time, make sure to check if it works properly with safety.

	Please be sure to read the instructions before use and follow as described.		
	Do not apply it to the person as following: 1. Who has acute disease or bloody disease over 38°C 2. Who has allergy and hemophilia 3. Who has melanosis of the face 4. Who has high internal heat 5. Who regularly visits hospital 6. Who takes medicines 7. Who is under 12 years old		
	Do not use it on the power other than AC 220V. It may cause fires or malfunction.		
	Do not use it near the place with high temperature, dust or humidity (Near the humidifier). It may cause malfunction.		
	Do not repair, disassemble or alter the controller and the body. It may cause malfunction.		
	Do not touch the device or detach the power cords with wet hands. It may cause electric shock.		
	Do not use the stripped cord or propped cord. It may cause fires(Be sure to replace it with another).		
	Do not use the propped cord which has crack or scares. It may cause scares on the skin.		
	When you don't use it for a long time, unplug it from the outlet.		
	Be sure to have an inspection for test. If any troubles are found, get it repaired.		
	No one other than the operator can touch it. It may cause accidents or malfunction.		
	Alternating current		Serial No
	Direct current		Manufacturer
	Date of manufacture		CE mark
	Attention, consult accompanying documents		Consult instruction for use

bio roller[®]
User Manual **G5**