

VELPEAU®



Thumb Wrist Brace VP0902C

[Product Name] : Thumb Wrist Brace (Drawstring Version)

[Product Code] :
 VP0902C-SL VP0902C-ML VP0902C-LL
 VP0902C-SR VP0902C-MR VP0902C-LR

[Product Specifications] :

Size	Product Code	Left/Right	Around Wrist(inch)
S	VP0902C-SL	Left	5"-6"
	VP0902C-SR	Right	
M	VP0902C-ML	Left	6"-8 1/2"
	VP0902C-MR	Right	
L	VP0902C-LL	Left	> 8 1/2"
	VP0902C-LR	Right	

[Product Features] :

- Anatomically shaped to provide functional support for the wrist while allowing free movement of the fingers.
- The functional drawstring closure method enhances the breathability and lightweight nature of the product. It allows for easy one-handed wearing and removal while providing the ability to quickly adjust the pressure on the wrist.
- Two aluminum alloy splints on the palm, along with elastic support strips on the back of the hand, provide effective control of palmar and dorsal flexion.
- The thumb splint is designed to be independent and removable, equipped with adjustable straps that secure the thumb in a functional position and restrict its movement. This allows for moderate flexibility while providing stability to the joint.

www.velpeau.com

- The thumb straps can be easily unwrapped at any time, providing a complete release for the thumb without the need to remove the entire brace. This feature ensures convenient care and maintenance.
- The brace is padded with sponge material for added comfort. It is also free from neoprene and latex, making it gentle on the skin and suitable for extended use.

[Product Componds] :

Polyester, Polyurethane, Cotton, Nylon, Polypropylene, Aluminum alloy, Polycarbonate.

[Application] :

- The wrist brace with thumb-spica splint can be used for all indications that require immobilisation of the wrist with thumb fixation, such as:
- Sprain, distorsion or stable fractures of the wrist and / or thumb.
 - Post-operative immobilization of the wrist and thumb.
 - De Quervain's Tenosynovitis.
 - Carpal tunnel syndrome (CTS).
 - Lesions of the collateral ligaments of the basal thumb joint.
 - Tendinitis and tenosynovitis of the wrist and / or thumb.
 - Radial nerve palsy of hand.
 - Inflammatory phase of rheumatic diseases.
 - After postoperative cast or splint removal.

www.velpeau.com

VELPEAU®



Shenzhen Weprotex Technology Product Co.,Ltd
 Add: 1710, Building 11, Phase II of Tian'an Yungu Industrial Park, Gangtuo Community, Bantian Street, Longgang District, Shenzhen, China
 Tel: 0755-89505533 E-mail: sales88@velpeau.com
 www.velpeau.com
 EC REP: No. 38, Boulevard Renard Benoit, 78680 Epone, France Made in China



[Instructions] :

- Proper sizing can be achieved by following the wear step diagram in the manual or watching the detailed wearing video on the official website www.velpeau.com.
- It is important to check the strap tightness to ensure stability. Avoid making it too tight, which can restrict blood flow and hinder wound healing, or loose, which may result in inadequate support.
- The product is designed for easy assembly and disassembly, requiring simple attachment or detachment of the hook and loop fastener.

[Precautions / Warnings] :

- Please consult a professional for guidance before using this product.
- This product is intended for adult and adolescent use only.
- It is designed for easy use, either independently or with the assistance of medical personnel or others when needed.
- If any discomfort occurs during use, adjust the tightness accordingly.
- The product is designed to be reusable, please wash the product with water before or after each use. Store it in a well-ventilated, dry place at room temperature.
- This product is designed for standalone use and should not be combined with other products. Its safety and effectiveness are not affected by external conditions.
- No specific special skills, knowledge, or training required for users. The product is made from materials that comply with the ROHS standard which is free from allergens or side effects on the human body.
- If there is an open wound on the skin, please consider bandaging the wound to prevent direct contact.

www.velpeau.com

[How to wear] :




www.velpeau.com

[Transportation And Storage Conditions] :

The product can be transported using general means of transportation and stored at room temperature. However, it must be fireproof, moisture-proof, dustproof, rodent-proof, and insect-proof.

[Washing] :

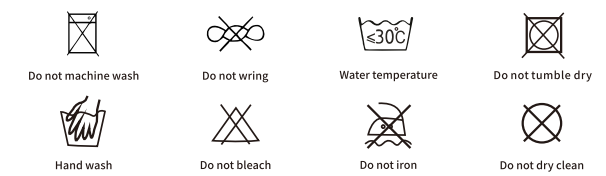

[Term of Validity] : 5 years

[Production Date] : See the package

[Manufacturer] : Shenzhen Weprotex Technology Product Co., Ltd.

[Production Address] : 1710, Building 11, Phase II of Tian'an Yungu Industrial Park, Gangtuo Community, Bantian Street, Longgang District, Shenzhen, China

[After-sale Service] : Shenzhen Weprotex Technology Product Co., Ltd.

[After-sale Service Address] : 1710, Building 11, Phase II of Tian'an Yungu Industrial Park, Gangtuo Community, Bantian Street, Longgang District, Shenzhen, China

[Postcode] : 518129

[Tel] : +86 - 755-89505533

[E-mail] : sales88@velpeau.com

[Website] : www.velpeau.com

[Date of Preparation] : Now.15,2022

[Modification Date] : Jul.15,2023

www.velpeau.com